


RELAXOR[®]

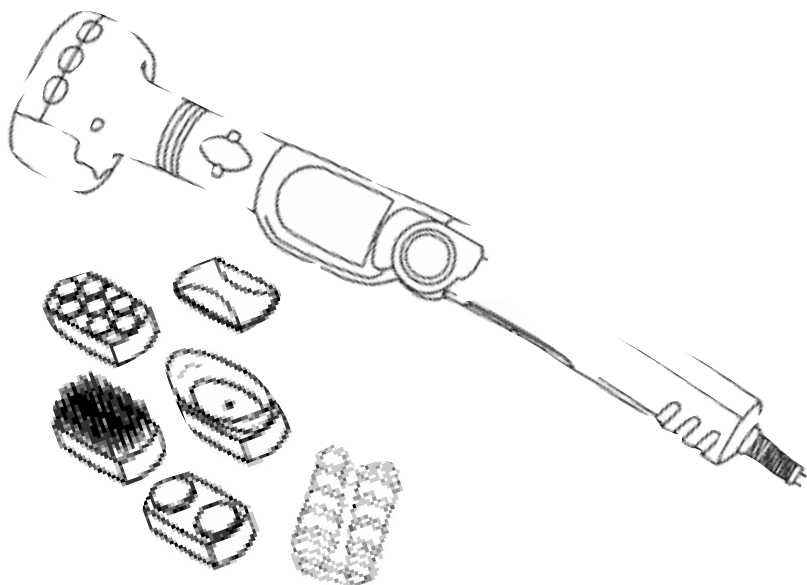

MASSAGER with HEAT

OWNER'S MANUAL &
OPERATING INSTRUCTIONS

Model MA717

REGENESIS[™]

IMPORTANT SAFEGUARDS

When using an electrical appliance, basic safety precautions should always be followed, including the following:

Read all instructions before using the appliance.

DANGER: To reduce the risk of electrical shock:

- Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- Do not reach for an appliance that has fallen into water. Unplug it immediately.
- Do not use while bathing or showering.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or any other liquid.
- Do not attempt to plug or unplug unit while you are standing in water.

WARNING: To reduce the risk of burns, fire, electric shock or injury to persons:

1. An appliance should NEVER be left unattended when plugged in. Unplug from electrical outlet when not in use.
2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.
3. The power cord must be fully extended during operation.
4. The Massager must never be immersed in water.
5. Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.
6. Use this appliance only for its intended use as described in this manual. Do not use Attachments not recommended by the manufacturer.
7. Do not stand on appliance. Never use while sleeping or falling asleep while using the Massager.
8. Ensure hands are dry when touching the Control Switch and plug.
9. Do not place or drop into water or any other liquid.
10. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped, damaged, or dropped in water. Please contact Consumer Service for return authorization.
11. Do not carry this appliance by power cord or use cord as a handle.
12. Keep cord away from heated surfaces.
13. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like. Never operate on a soft surface such as a bed or couch where openings may be blocked.
14. Never drop or insert any object into any opening.
15. Do not use outdoors.
16. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

IMPORTANT SAFEGUARDS (Con't.)

17. To disconnect, push all Controls so that they are in the OFF position, then remove plug from outlet.
18. This appliance is designed for household use only.

CAUTION: Please read all instructions carefully before operating.


- If you have any concerns regarding your health, consult your doctor before using the Hand Massager.
- Individuals with pacemakers should consult a physician before using this product.
- Never leave the appliance unattended, especially if children are present.
- Never cover the appliance when it is in operation.
- Do not use this product for more than 20 minutes at a time.
- Extended use for longer periods of time could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- Never use this appliance directly on swollen or inflamed areas or skin eruptions.
- Never use this product in bed.
- This appliance should never be used by an individual suffering from any physical ailment that would limit the user's capacity to operate the Controls.
- This unit should not be used by children without adult supervision.
- In case of pregnancy, diabetes or illness, consult your doctor before using the Massager.
- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.
- Never use any Massager on open wounds, discolored areas, or any area of the body that is swollen, burned, inflamed or where skin eruptions or sores are present.
- Do not use on legs with varicose veins.
- Massage should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.

SAVE THESE INSTRUCTIONS FOR HOUSEHOLD USE ONLY

POLARIZED PLUG

This appliance has a **polarized plug** (one blade is wider than the other). To reduce the risk of electric shock, this plug is intended to fit into a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.

Congratulations on your purchase of your RegenesiTM Massager with Heat, one of the latest in a long line of personal care products brought to you by Relaxor.[®] Like our other fine products, the Regenesi Massager is engineered for quality, comfort and safety. Its powerful motor, 6 Massage Heads and Heat option combine to give you a wonderfully soothing experience. It may provide penetrating relief and relaxation for sore and tense muscles. This massaging product will also serve as a great addition to your health and fitness program.

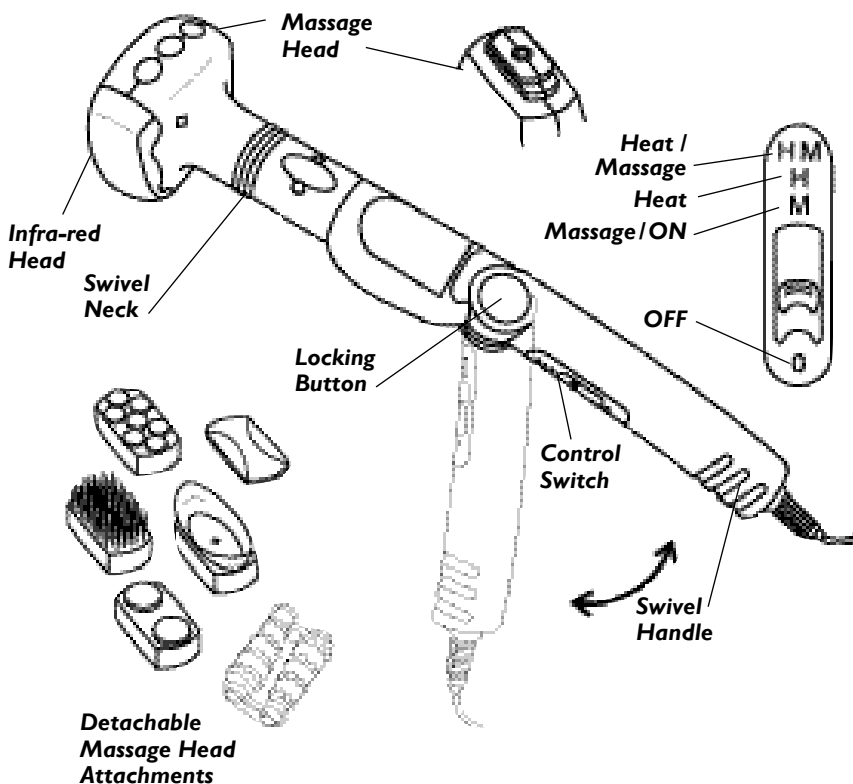

In order to enhance your enjoyment and maximize the benefits of your purchase, please take a few moments to read the following instructions.

INSTRUCTIONS FOR USE

1. Plug unit into AC electrical outlet. Make sure unit is set to "OFF" position before inserting plug into electrical outlet.
2. To turn unit ON, slide Control Switch up one position to the **MASSAGE/ON** position indicated with an "M." To use **HEAT ONLY**, slide the Control Switch up one position further to the **HEAT** position indicated with an "H." To use **HEAT and MASSAGE** together, slide the Control Switch all the way up to the top **HEAT/MASSAGE** position indicated by the "HM."
3. Use the Massager on arms, legs, back or any area you wish to treat. Slowly move the unit around for maximum comfort.
4. When finished, slide the Control Switch to the "OFF" position and unplug the unit.

REGENESIS™ MASSAGER ATTACHMENTS



**Massage
Rollers**
(P/N 50542)



Suction
(P/N 50543)



**Double
Action**
(P/N 50544)



Scalp
(P/N 50545)



Multi-Point
(P/N 50546)



Deep Muscle
(P/N 50547)

Massage Rollers

For all-over body relaxation, additional stimulation is provided by the dimpled rollers

Suction

Makes use of a slight vacuum. May be beneficial for skin, thighs and buttocks.

Double Action

For thighs, calves, shoulders and large muscle areas

Scalp

To stimulate for healthy hair, or use at the nape of the neck to help ease tension

Multi-Point

For general purpose use, and all-over massage

Deep Muscle

For specific parts of muscles and small areas such as hands and feet.

TO CHANGE ATTACHMENTS

Your Massager comes equipped with 6 Custom Massage Head Attachments.

- 1. Unplug Massager from outlet before putting on or taking off Attachments.*
- 2. To attach, simply push the Attachment into place on the Massage Head.*

SWIVEL NECK

Gripping the Handle firmly, twist the Neck to bring either the Heat Lamp side or Massage Head side into the position required. The Neck does not turn completely around, only 270 degrees. It has 4 pre-set positions, 90 degrees apart.

ADJUSTABLE HANDLE

The Regensis™ Massager's Handle bends to 4 pre-set positions to help massage those awkward areas. Gripping the upper part of the Massager firmly, press the Locking Button and push the Handle into the position required. Release the Locking Button when it is at the desired location and it will click into place.

MAINTENANCE AND CLEANING

Unplug the unit and allow it to cool. Use a soft, slightly damp cloth to wipe. Never allow water or any liquids to come into contact with the handle. Do not use any chemicals or abrasive cleaners.

Periodically untangle cord if it becomes twisted.

This unit has no user serviceable parts. If unit is damaged or malfunctions, refer to the warranty section of this manual.

STORAGE

Store unit in its box or in a clean, dry place.

Do not wrap cord around unit as this may cause it to break.

LIMITED ONE YEAR WARRANTY

Warranty: This Relaxor® product is warranted by Salton, Inc. to be free from defects in materials or workmanship for a period of (1) year from the original purchase date. This product warranty covers only the original consumer purchaser of the product.

Warranty Coverage: This warranty is void if the product has been damaged by accident in shipment, unreasonable use, misuse, neglect, improper service, commercial use, repairs by unauthorized personnel, normal wear and tear, improper assembly, installation or maintenance abuse or other causes not arising out of defects in materials or workmanship. This warranty is effective only if the product is purchased and operated in the USA, and does not extend to any units which have been used in violation of written instructions furnished by with the product or to units which have been altered or modified or to damage to products or parts thereof which have had the serial number removed, altered, defaced or rendered illegible.

Implied Warranties: ANY IMPLIED WARRANTIES WHICH THE PURCHASER MAY HAVE ARE LIMITED IN DURATION TO ONE (1) YEAR FROM THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

Warranty Performance: During the above one-year warranty period, a product with a defect will be either repaired or replaced with a reconditioned comparable model (at our option) when the product is returned to the Repair Center, or the purchase price refunded. The repaired or replacement product will be in warranty for the balance of the one-year warranty period and an additional one-month period. No charge will be made for such repair or replacement.

Service and Repair: Should the appliance malfunction, you should first call toll-free 1 (800) 233-9054 between the hours of 9:00 am and 5:00 pm Central Standard Time and ask for CONSUMER SERVICE stating that you are a consumer with a problem. Please refer to model number MA717 when you call.

In-Warranty Service (USA): For an appliance covered under the warranty period, no charge is made for service or postage. Call for return authorization 1 (800) 233-9054.

Out-of Warranty Service: A flat rate charge by model is made for out-of-warranty service. Include \$6.00 (U.S.) for return shipping and handling. We will notify you by mail of the amount of the charge for service and require you to pay in advance for the repair or replacement.

For Products Purchased in the USA, but Used in Canada: You may return the product insured, packaged with sufficient protection, and postage and insurance prepaid to the USA address listed below. Please note that all customs duty / brokerage fees, if any, must be paid by you and we will require you to pay the cost of customs duty / brokerage fees to us in advance of our performing any service.

Risk During Shipment: We cannot assume responsibility for loss or damage during incoming shipment.

For your protection, carefully package the product for shipment and insure it with the carrier. Be sure to enclose the following items with your appliance: any accessories related to your problem, your full return address and daytime phone number, a note describing the problem you experienced, a copy of your sales receipt or other proof of purchase to determine warranty status. C.O.D. shipments cannot be accepted.

To return the appliance, ship to:

ATTN: Repair Center
Salton, Inc.
708 South Missouri Street
Macon, MO 63552

To contact us, please write to, call, or email:

Salton, Inc.
P.O. Box 1637
Columbia, MO 65205-1637
1-800-233-9054
E-mail: Salton@Saltonusa.com

Limitation of Remedies: No representative or person is authorized to assume for Salton, Inc. any other liability in connection with the sale of our products. There shall be no claims for defects or failure of performance or product failure under any theory of tort, contract or commercial law including, but not limited to negligence, gross negligence, strict liability, breach of warranty and breach of contract. Repair, replacement or refund shall be the sole remedy of the purchaser under this warranty, and in no event shall Salton be liable for any incidental or consequential damages, losses or expenses. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

Legal Rights: This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

For more information on Salton products: visit our website: www.eSalton.com

©2003 Salton, Inc. All rights reserved.

Made in China

P/N 61276

2/03