

RELAXOR®



DEEP RELIEF™ 9 MASSAGE MOTORS & HEAT

OWNER'S MANUAL &
OPERATING INSTRUCTIONS

Model FBM969P

FULL BODY MAT

IMPORTANT SAFEGUARDS

When using your RELAXOR® SEAT TOPPER, basic safety precautions should always be followed, including the following:

1. Read all instructions before using.

To reduce risk of electric shock:

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not reach for a product that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.

WARNING! OPERATE YOUR SYSTEM ACCORDING TO THIS MANUAL. USE ONLY THE POWER DEVICE SUPPLIED

To reduce risk of burns, fire, electric shock, or injury to persons:

1. Your RELAXOR® product should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate the RELAXOR® product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact your dealer for instructions on returning it for examination and repair.
4. Close supervision is necessary when this product is used by, on, or near children or disabled persons.
5. Do not use outdoors.
6. Keep the cord away from heated surfaces.
7. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
8. Use heated surfaces carefully. May cause serious burns. Do not use heated surfaces over insensitive skin areas or in the presence of poor circulation.
9. If you are under a physician's care and/or have a condition requiring the use of any medical device, please consult your health care provider prior to using this product.


SAVE THESE INSTRUCTIONS FOR HOUSEHOLD USE ONLY

OPERATING INSTRUCTIONS

1. Control Wand is permanently attached.
2. Leave the Power Switch on the Control Wand in the OFF position. Put the RELAX/BOTH/HEAT Switch in the RELAX position.
3. Plug the cord from the Transformer into the remaining receptacle coming from your product, and then plug the Transformer Body into the electrical wall outlet. (120V AC - 60 Hz power.)
4. Sit on your Massaging Mat and familiarize yourself with the system operation by following the next instructions.
5. Place the WAVE/ALL/ZONE Switch in the ZONE position. Slide the INTENSITY Switch to LO. Turn on the Power Switch. You can now manually position the massage action to any one of five body locations. Each time you press the SELECT Button you will move the massage action one position up or down, as indicated by the green LED light on the Control Wand. Choose the desired intensity by selecting either LO, MED, or HI.
6. Set WAVE/ALL/ZONE Switch to WAVE. Slowly slide the SPEED Knob up and down to find the wave action you like best.

The green lights on the Control Wand will flash in sequential order from top to bottom, then back up to top continuously repeating this cycle. The lights show the activation of the corresponding pulsing transducers and you will feel a gentle massaging action at the positions shown by the lights.

7. Move the WAVE/ALL/ZONE Switch to ALL for Massaging Motors to alternate between zones 1 and 2.
8. Move the RELAX/BOTH/HEAT Switch to HEAT, and after a short warmup time, enjoy comforting warmth in the lower back (lumbar) area. Now move the Switch to BOTH and indulge in the ultimate combination of Relaxor® massage and heat.
9. Select any combination of the above described modes, intensities and heat.

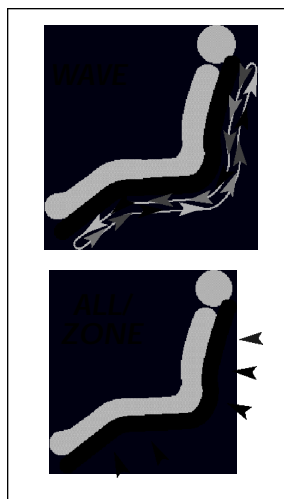

OPERATING INSTRUCTIONS (Continued)

You have now completed all of the operational choices provided on the Control Wand.

The system will automatically shut-down after approximately 15-minutes of operation. To restart, simply move the POWER Switch to OFF then back to ON.

Now close your eyes and enjoy your Relaxor® Deep Relief™ Full Body Mat.

NOTE: It is normal for the Control Wand to become warm during use.

If your Relaxor® System does not operate properly, or you need assistance for any reason, call Consumer Service.

USER MAINTENANCE INSTRUCTIONS

No user maintenance of this product is required. If the System does not function, be sure...

(a) The Transformer is plugged into an electrical outlet.

(b) The Transformer cord is plugged into, and fully engaged, in the receptacles.


If there is a malfunction of any kind, contact your dealer for instructions on returning it for examination and repair.

CLEANING

You may clean the massaging system with a damp cloth and warm water.

DO NOT SATURATE. Allow adequate time to dry.

CONTROL WAND INSTRUCTIONS

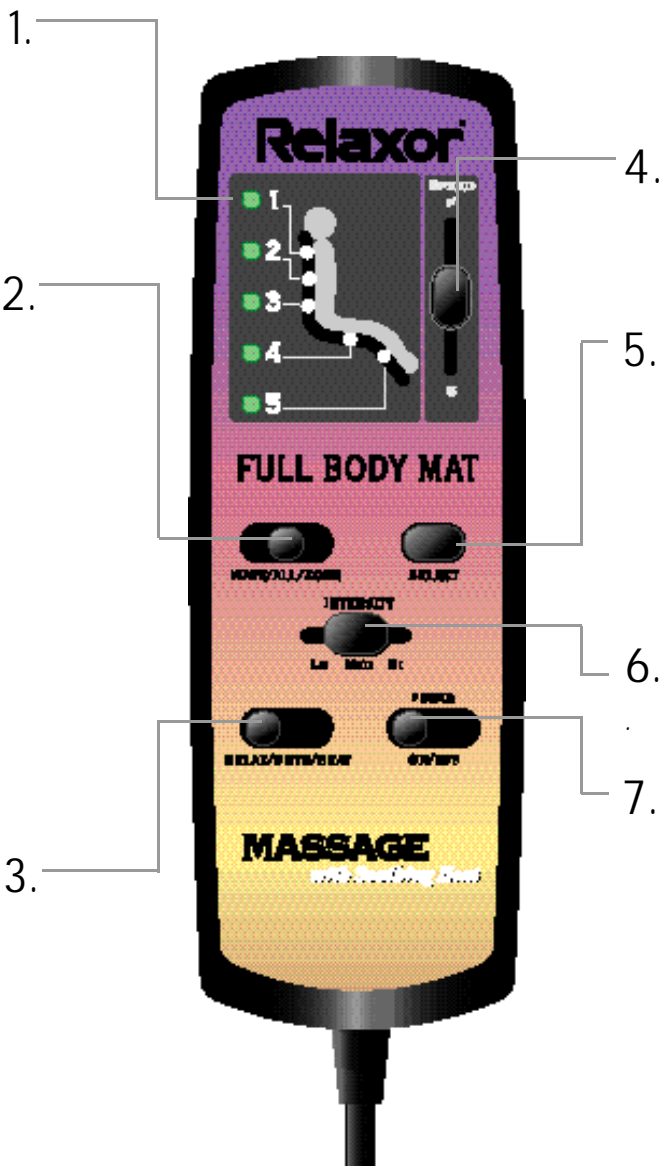

CONTROL WAND INSTRUCTIONS (Continued)

1.

ZONE 1: *Upper Shoulders*

ZONE 2: *Mid back area*

ZONE 3: *Lumbar area*

ZONE 4: *Thighs*

ZONE 5: *Calves*

2.

WAVE/ALL/ZONE:

WAVE - For a wave-like sensation slide the WAVE/ALL/ZONE Button to WAVE. (Adjust intensity and speed to your preference.)

ALL - For an overall massage slide the WAVE/ALL/ZONE Button to ALL.

ZONE - To isolate individual massaging zone, slide the WAVE/ALL/ZONE Button to ZONE (Adjust intensity and speed to your preference.)

3.

RELAX/BOTH/HEAT:

RELAX: For massage without heat slide the RELAX/BOTH/HEAT Button to RELAX.

BOTH: For massage with heat slide the RELAX/BOTH/HEAT Button to BOTH.

HEAT: For heat without massage slide the RELAX/BOTH/HEAT Button to HEAT.

4.

SPEED: Adjust from down to up to control the cycling speed of the Massaging Motors. Lighted LEDs will show what region of the body is active with massage.

5.

SELECT: Press the SELECT Button to manually position the massage action to any one of five body locations. The massage action moves one position up or down, as indicated by the green LED LIGHT.

6.

INTENSITY: Choose LO/MED/HI for desired intensity.

7.

POWER:

- Turn on: Slide OFF/ON Button to ON.
- Turn off: Slide OFF/ON Button to OFF.

15 MINUTE TIMER: To reset after Timer shuts off, slide the POWER Button to OFF, then back to ON.

(NOTE: The green LED light will show the zone that is active).

LIMITED ONE YEAR WARRANTY

Warranty: This Relaxor® product is warranted by Salton, Inc. to be free from defects in materials or workmanship for a period of (1) year from the original purchase date. This product warranty covers only the original consumer purchaser of the product.

Warranty Coverage: This warranty is void if the product has been damaged by accident in shipment, unreasonable use, misuse, neglect, improper service, commercial use, repairs by unauthorized personnel, normal wear and tear, improper assembly, installation or maintenance abuse or other causes not arising out of defects in materials or workmanship. This warranty is effective only if the product is purchased and operated in the USA, and does not extend to any units which have been used in violation of written instructions furnished by with the product or to units which have been altered or modified or to damage to products or parts thereof which have had the serial number removed, altered, defaced or rendered illegible.

Implied Warranties: ANY IMPLIED WARRANTIES WHICH THE PURCHASER MAY HAVE ARE LIMITED IN DURATION TO ONE (1) YEAR FROM THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

Warranty Performance: During the above one-year warranty period, a product with a defect will be either repaired or replaced with a reconditioned comparable model (at our option) when the product is returned to the Repair Center, or the purchase price refunded. The repaired or replacement product will be in warranty for the balance of the one-year warranty period and an additional one-month period. No charge will be made for such repair or replacement.

Service and Repair: Should the appliance malfunction, you should first call toll-free 1 (800) 233-9054 between the hours of 9:00 am and 5:00 pm Central Standard Time and ask for CONSUMER SERVICE stating that you are a consumer with a problem. Please refer to model number FBM969P when you call.

In-Warranty Service (USA): For an appliance covered under the warranty period, no charge is made for service or postage. Call for return authorization 1 (800) 233-9054.

Out-of-Warranty Service: A flat rate charge by model is made for out-of-warranty service. Include \$8.00 (U.S.) for return shipping and handling. We will notify you by mail of the amount of the charge for service and require you to pay in advance for the repair or replacement.

For Products Purchased in the USA, but Used in Canada: You may return the product insured, packaged with sufficient protection, and postage and insurance prepaid to the USA address listed below. Please note that all customs duty / brokerage fees, if any, must be paid by you and we will require you to pay the cost of customs duty / brokerage fees to us in advance of our performing any service.

Risk During Shipment: We cannot assume responsibility for loss or damage during incoming shipment. For your protection, carefully package the product for shipment and insure it with the carrier. Be sure to enclose the following items with your appliance: any accessories related to your problem, your full return address and daytime phone number, a note describing the problem you experienced, a copy of your sales receipt or other proof of purchase to determine warranty status. C.O.D. shipments cannot be accepted.

To return the appliance, ship to:

ATTN: Repair Center
Salton, Inc.
708 South Missouri Street
Macon, MO 63552

To contact us, please write to, call, or email:

Salton, Inc.
P.O. Box 1637
Columbia, MO 65205-1637
1-800-233-9054
E-mail: Salton@Saltonusa.com

Limitation of Remedies: No representative or person is authorized to assume for Salton, Inc. any other liability in connection with the sale of our products. There shall be no claims for defects or failure of performance or product failure under any theory of tort, contract or commercial law including, but not limited to negligence, gross negligence, strict liability, breach of warranty and breach of contract. Repair, replacement or refund shall be the sole remedy of the purchaser under this warranty, and **in no event shall Salton be liable for any incidental or consequential damages, losses or expenses.** Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

Legal Rights: This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

For more information on Salton products: visit our website: www.salton-maxim.com

©2001 Salton, Inc. All rights reserved.

Made in China

P/N 60753

06/01