


Relaxor[®]



(((9)))Massaging Motor Full Body Mat with Flip Pillow and Soothing Heat

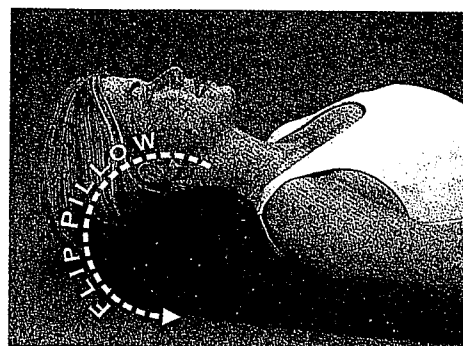

Model
989

OWNERS MANUAL and OPERATING INSTRUCTIONS

Dear Customer:

Thank you for purchasing a **Relaxor** Full Body Mat with Flip Pillow. This unique massage product is a scientifically designed electro-mechanical system with built-in flip pillow and multiple pulsing transducers. It provides many customized settings allowing you to create your own personalized massage.

Total relaxation through finger-like pressure massage is now yours to command. With built-in flip pillow and soothing heat, you have a winning combination for your comfort.

In order to enhance your enjoyment and maximize the benefits of your new purchase, please take a few moments to read the following instructions.

Relaxor[®]

FEATURES AND BENEFITS




Experience Deep Penetrating Relief

ENJOY THESE PRODUCT BENEFITS!

Full Body Mat with Flip Pillow and Soothing Heat

- Dozens of programmable features.
- Great for recliner, bed or sofa.
- Membrane microprocessor controlled wand.
- Fashionable design for ease of control.
- 15/30 minute shut-off timer.

- Built in adjustable flip pillow cradles the neck offering optimum support and comfort (can be used with or without).
- Durable high quality foam.
- Low voltage motors enhance safety.
- Comes with AC Adaptor.

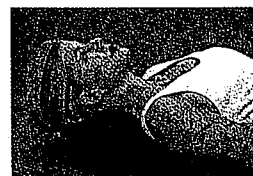


Relaxor®

How To Operate Your Full Body Mat with Flip Pillow

POWER:

- Turn on: Press PWR button
- Turn off: Press PWR button again
- 15 and 30 minute timers:**
- Press PWR button:
Green light = 15 minute timer
- Hold down PWR button:
Amber light = 30 minute timer
- To reset after timer shuts off: Press the PWR button
- To reset timer while on: Press PWR button twice

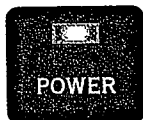

SOOTHING HEAT:


LO or HI

- Works independently of massage
- 30 minute timer
- To reset after timer shuts off:
Press the LO or HI button
- To reset timer while on:
turn HEAT off, then turn HEAT back on



((9)) Motors in 5 strategic zones

ZONES:

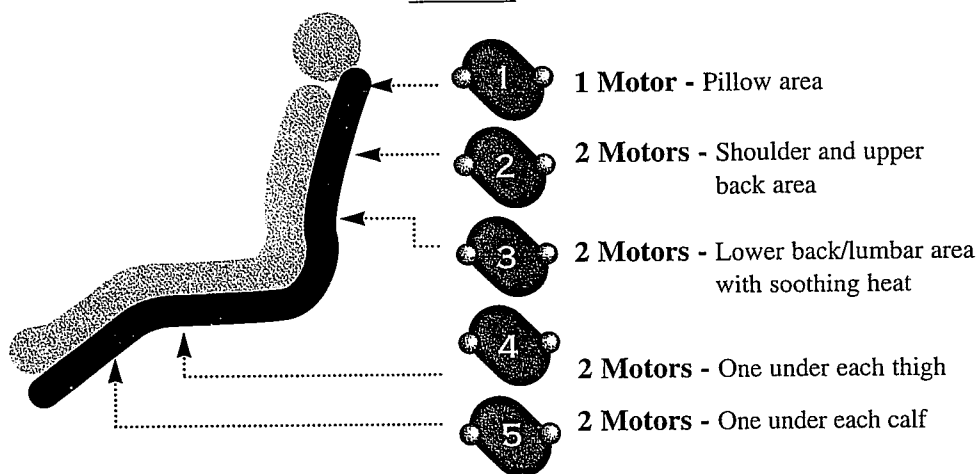


7 Different Styles of Massage

Note: Speed and Intensity are regulated by the SPEED and INTENSITY buttons (refer to SPEED and INTENSITY)

Massaging motors cycle down and up in a wave-like sensation

- Press WAVE
- Select any or all Zones (1-5)
(To deactivate zone, re-press Zone Key (1-5))

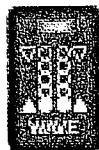


Massaging motors pulsate on and off

- Press PULSE
- Select any or all Zones (1-5)
(To deactivate zone, re-press Zone Key (1-5))



Isolates massage area

- Press SELECT
- Select any or all Zones (1-5)
(To deactivate zone, re-press Zone Key (1-5))

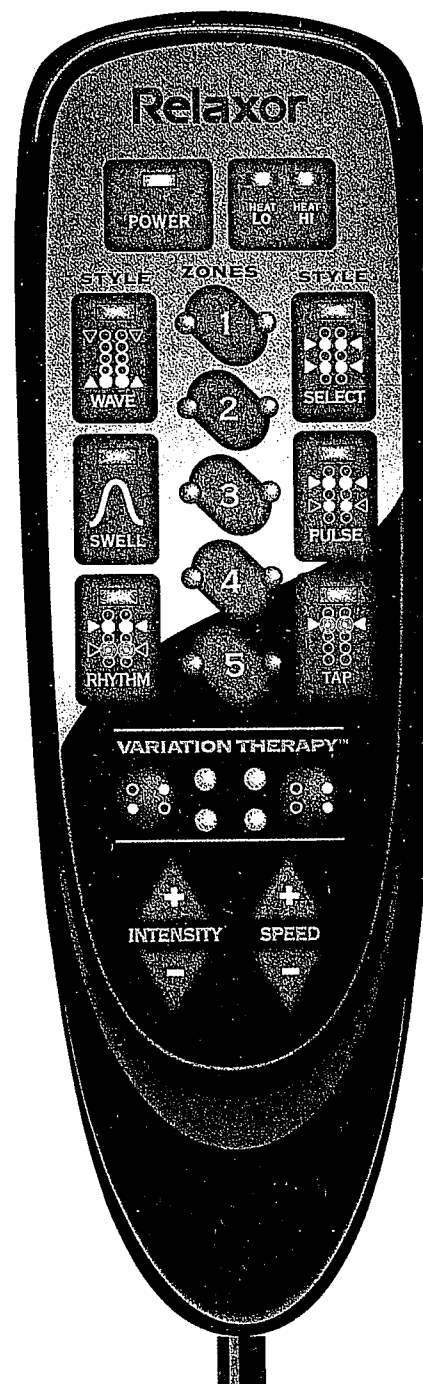

WAVE

PULSE


SELECT

Relaxor®

How To Operate Your Full Body Mat with Flip Pillow



SWELL

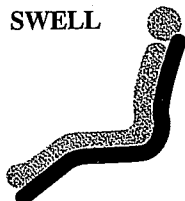

Massaging motors cycle up and down in intensity

- Press SWELL
- Select any or all zones (1-5)
(To deactivate zone, re-press Zone Key (1-5))
- To adjust the peak intensity of the SWELL feature, use the INTENSITY + / - buttons
- To adjust the cycle time of the SWELL feature, use the SPEED + / - buttons



TAPPING



Massaging motors are in tapping mode

- Press TAP
- Select any or all zones (1-5)
(To deactivate zone, re-press Zone Key (1-5))
- To adjust the speed of the TAPPING, use the SPEED + / - buttons
- To adjust the intensity of the TAPPING, use the INTENSITY + / - buttons



RHYTHM



Massaging motors will alternate in tapping and select continuous mode

- Press RHYTHM
- Select any or all zones (1-5)
(To deactivate zone, re-press Zone Key (1-5))
- To change the duration of the RHYTHM Mode, use the SPEED + / - buttons
- To adjust the intensity of the RHYTHM, use the INTENSITY + / - buttons



VARIATION THERAPY:

Two pre-set variations using the top four massaging motors in Zones 2 & 3.
(Located in shoulder/upper back and lower back areas)



Note: Speed and Intensity are regulated by the SPEED and INTENSITY buttons (refer to SPEED and INTENSITY)

MASSAGE MOTOR ACTION:

Variation Therapy lights show activated motors

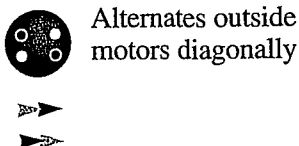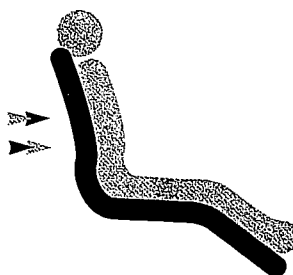

Alternates outside motors diagonally

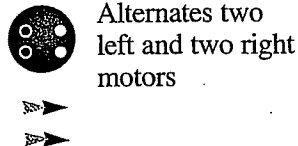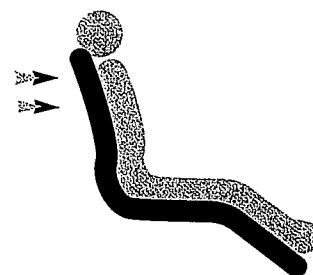

Alternates two left and two right motors

INTENSITY:

- Controls strength of massage
- Press + / - keys for desired effect



SPEED:

- Controls cycle speed of massaging motors in WAVE, PULSE, SWELL, RHYTHM, TAP, or VARIATION THERAPY
- SPEED button does not apply in SELECT (continuous) mode
- Press + / - keys for desired effect



IMPORTANT SAFETY INSTRUCTIONS

Read all instructions before using the **Relaxor** Full body Mat with Flip Pillow

DANGER!

To reduce risk of electric shock:

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not reach for a product that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.

WARNING!

To reduce risk of burns, fire, electric shock, or injury to persons:

1. Your Full Body Mat product should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate the Full Body Mat product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact your dealer for instructions on returning it for examination and repair.
4. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
5. Do not use outdoors.
6. Keep the cord away from heated surfaces.
7. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
8. Use heated surfaces carefully. May cause serious burns. Do not use heated surfaces over insensitive skin areas or in the presence of poor circulation.
9. If you are under a physician's care and/or have a condition requiring the use of any medical device, please consult your health care provider prior to using this product.

MAINTENANCE AND CARE INSTRUCTIONS:

No user maintenance of this product is required. If the System does not function, be sure...

- (a) The TRANSFORMER is plugged into an electrical outlet.
- (b) The Transformer cord and COMMAND WAND cord are plugged into, and fully engaged, in their respective receptacles.

If there is a malfunction of any kind, contact your dealer for instructions on returning it for examination and repair.

How To Clean:

You may clean the seat topper with a damp cloth and warm water.

DO NOT SATURATE. Allow adequate time to dry.


SAVE THESE INSTRUCTIONS!

OPERATE
YOUR
SYSTEM
ACCORDING
TO THIS
MANUAL


USE ONLY
THE POWER
DEVICE
SUPPLIED

Relaxor®


OTHER FINE PRODUCTS FROM THE LEADER IN RELAXATION

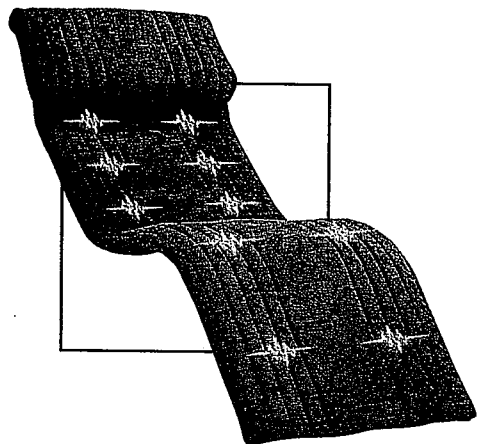

#824 CONTOURCARESS™
10 soft tissue massaging motors plus
soothing heat. Contour pillow built in.


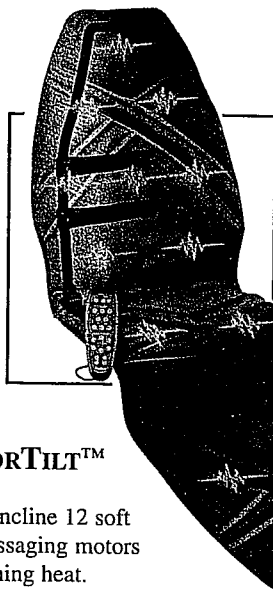


**#998
COMFORTILT™
PLUS**
Variable incline 12 soft
tissue massaging motors
plus soothing heat.
Use with or without frame.

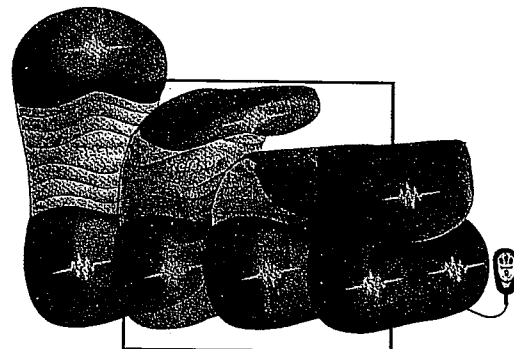

#340 BACK2SHOULDERS™
Converts in seconds from back massager
to neck and shoulder massager. 3 motors,
18 customized settings and soothing
heat. Great for home, office and auto.

**#721
SOUND & LIGHT THERAPY
WITH ALARM CLOCK**
Six soothing environmental sounds
to promote rest and relaxation. Full-
featured clock, headphone outlet and
comforting soft illuminating light.
Sound & Light Therapy performs as
a friendly wake up call.

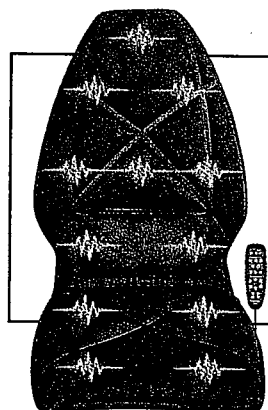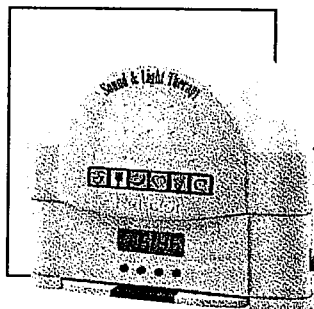

**#964 RELAXOR®
MASSAGING SEAT TOPPER**
12 Motor Massage Seat Topper with
Soothing Heat. Features dozens of
settings for an all over massage.
Great for home, office and auto.

Relaxor®

The Leader in Relaxation™

Relax and Enjoy!

Limited One Year Warranty

If your FULL BODY MAT WITH FLIP PILLOW does
not operate properly, or you need assistance for
any reason, call:

Customer Service - JB Research, Inc.
Los Angeles Corporate Office • Bellflower, CA
1-800-RELAXOR

Outside U.S. write to:
Customer Service
9062 Rosecrans Ave. Bellflower, CA 90706
© Copyright 1998 JB Research, Inc.
PATENTS PENDING



(Adaptor)

#15155

899303-00