


Dr. Seuss'
**HOW THE
GRINCH**
STOLE
CHRISTMAS!

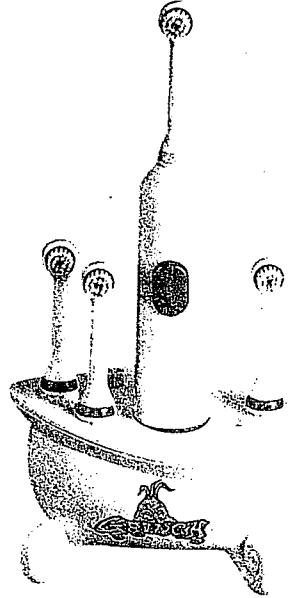

**WHOVILLE
RECHARGEABLE TOOTHBRUSH**

MODEL: GRIN8TB

USE & CARE GUIDE

**salton®**

IMPORTANT SAFEGUARDS

When using electrical appliances, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER - To reduce the risk of electrocution:

1. Do not use while bathing.
2. Do not place or store product where it can fall or be pulled into a tub or sink.
3. Do not place in or drop into water or other liquid.
4. Do not reach for a product that has fallen into water. Unplug it immediately.

WARNING: To reduce the risk of burns, electrocution, fire, or injury to persons:

1. A product should never be left unattended when plugged in.
2. Close supervision is necessary when this appliance is used by, on or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged in any manner or dropped into water. Return the product to a Service Center for examination and repair.
5. Keep the cord away from hot or heated surfaces.
6. Never use while drowsy.
7. Never drop or insert any object into any opening.
8. Do not use outdoors or operate where aerosol (spray) products are being used. Do not operate where oxygen is being administered.
9. The supply cord of this appliance cannot be replaced; if the cord is damaged, the appliance should be discarded.

THIS APPLIANCE IS FOR HOUSEHOLD USE ONLY
SAVE THESE INSTRUCTIONS

POLARIZED PLUG

This appliance has a polarized plug (one blade is wider than the other). To reduce the risk of electric shock, this plug is intended to fit into a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.

SHORT CORD INSTRUCTIONS


A short power-supply cord is provided to reduce the risk resulting from becoming entangled in or tripping over a longer cord.

Longer extension cords are available and may be used if care is exercised in their use.

If an extension cord is used, (1) the marked electrical rating of the extension cord should be at least as great as the electrical rating of the appliance, and (2) the longer cord should be arranged so that it will not drape over the counter top or table top where it can be pulled on by children or tripped over unintentionally.

OPERATING INSTRUCTIONS

PRODUCT MAY VARY SLIGHTLY FROM ILLUSTRATIONS

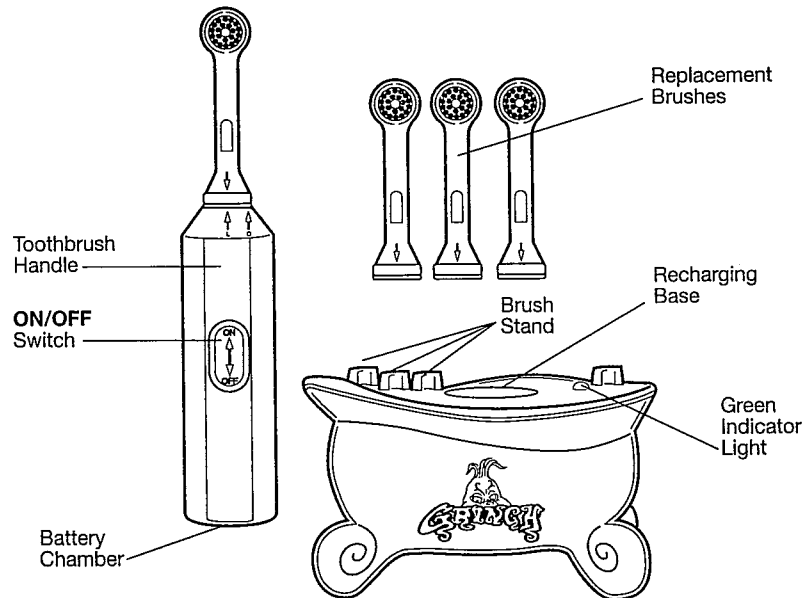

HOW TO CHARGE YOUR TOOTHBRUSH


Your Grinch Toothbrush requires charging for 12 hours before use. This product can only be used in conjunction with the Recharging Base supplied.

1. Make sure the Power Switch is in the **OFF** position.
2. Place the Toothbrush into the Recharging Base.
3. Plug the cord on the Recharging Base into a 120V electrical outlet.
4. Charging is indicated by the illumination of the Green Indicator Light on the Recharging Base. The normal charging time for this product is 12 hours.

The fully charged Toothbrush will have approximately 2 hours of operating time before recharging is necessary.

INSTRUCTIONS FOR USE

After the Toothbrush has been charged for 12 hours it is ready for use. Brushes are interchangeable and color coded for family use.

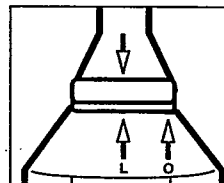


1. Attach one of the four color coded Brushes onto the Toothbrush Handle by grasping the Brush by the neck, not the bristles, and placing the Brush onto the top of the Toothbrush Handle. Align the arrow on the Brush to the "O" (Open) arrow on the Handle. Twist the Brush to the left until it aligns with the "L" (Lock) arrow on the Handle.
2. Place toothpaste onto the Brush. It is recommended that you consult your dentist as to which type of toothpaste is best for you.
3. Position yourself over the sink and direct the Toothbrush against your teeth. Slide the Switch to **ON** and begin brushing. Please consult your dentist for the appropriate brushing pattern to follow for your dental care.
4. When brushing is complete, slide the Switch to **OFF** and return the Toothbrush to the Recharging Base.
5. Remove the Brush by twisting it to the "O" (Open) arrow on the Toothbrush Handle. Rinse Brush off under running water and replace onto Brush Stand.
6. Plug the cord on the Recharging Base into a 120V electrical outlet if recharging is needed.
7. When not recharging the unit, always unplug the cord from the outlet.

REPLACEMENT BRUSHES


It is recommended to change the Brushes every 3 months to maintain cleaning efficiency and proper dental health. Replacement Brushes (in packs of two) can be purchased through our Consumer Service Department by calling 1-800-233-5054 or 1-888-889-0899. Please refer to Model Number GRIN8RTB when you call.

HOW TO REPLACE THE BATTERIES

The Battery Chamber is located on the bottom of the Toothbrush. Please refer to the diagram below for instructions on how to remove the Battery Chamber Cover. Insert 2 "AA" rechargeable batteries (replacement batteries not included) as shown on the diagram inside the Battery Chamber. Replace the Battery Chamber Cover.

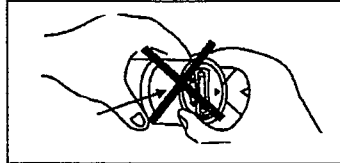

Do not pull or yank on the Toothbrush Handle Base.

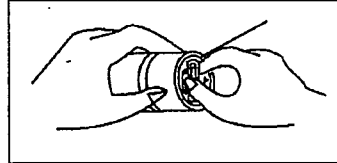

For easy battery removal, turn the Battery Chamber Cover Knob from "L" (LOCK) to "O" (OPEN).

USER MAINTENANCE INSTRUCTIONS

This appliance contains no user serviceable parts. Do not try to repair it yourself. Refer service, if necessary, to qualified service personnel.

TO CLEAN: Unplug appliance and allow it to cool. Use a soft, slightly damp cloth or sponge to wipe the surface. Never use abrasive cleaners or solvents, as they may scratch the surface.

TO STORE: Unplug appliance and allow it to cool. Store unit in its box or in a clean, dry place. Never store unit while it is still plugged in. Do not put any stress on the cord where it enters the unit, as it could cause the cord to fray and break.

ONE-YEAR LIMITED WARRANTY

This Salton, Inc. product warranty extends to the original consumer purchaser of the product.

Warranty Duration: This product is warranted to the original consumer purchaser for a period of one (1) year from the original purchase date.

Warranty Coverage: This product is warranted against defective materials or workmanship. This warranty is void if the product has been damaged by accident, in shipment, unreasonable use, misuse, neglect, improper service, commercial use, repairs by unauthorized personnel or other causes not arising out of defects in materials or workmanship. This warranty is effective only if the product is purchased and operated in the USA, and does not extend to any units which have been used in violation of written instructions furnished by Salton, Inc., or to units which have been altered or modified or to damage to products or parts thereof which have had the serial number removed, altered, defaced or rendered illegible.

Warranty Disclaimers: This warranty is in lieu of all warranties expressed or implied and no representative or person is authorized to assume for Salton, Inc. any other liability in connection with the sale of our products. There shall be no claims for defects or failure of performance or product failure under any theory of tort, contract or commercial law including, but not limited to, negligence, gross negligence, strict liability, breach of warranty and breach of contract. Some states do not allow the exclusion or limitations of implied warranties or consequential damages, so the above limitations or exclusions may not apply to you.

Salton, Inc. is not responsible or liable for indirect, special or consequential damages arising out of or in connections with the use or performance of the product or other damages with respect to loss of property, or loss of revenues or profit.

Legal Remedies: This warranty gives you specific legal rights, and you may also have other rights that vary from state to state.

Warranty Performance: During the above one-year warranty period, a product with a defect will be either repaired or replaced with a reconditioned comparable unit (at Salton, Inc.'s option) when the product is returned to the Salton, Inc. facility listed below. The repaired or replacement product will be in warranty for the balance of the one-year warranty period and an additional one-month period. No charge will be made for such repair or replacement.

Service and Repair

If service is required, you should first call toll-free 1-800-233-9054 **between the hours of 8:00 am and 5:00 pm Central Standard Time** and ask for **CONSUMER SERVICE**. Please refer to Model GRIN8TB when you call.

In-Warranty Service for a product covered under the warranty period, no charge is made for service and return postage. Please return the product insured, packaged with sufficient protection, and postage and insurance prepaid to the address listed below. Call for return authorization.

Out-of-Warranty Service: There will be charges rendered for repairs made to the product after the expiration of the aforesaid one (1) year warranty period, after purchaser is advised appropriately. Include \$6.00 (U.S.) for return shipping and handling.

Salton, Inc. cannot assume responsibility for loss or damage during incoming shipment. For your protection, carefully package the product for shipment and insure it with the carrier. Be sure to enclose the following items with your unit: any accessories related to your problem; your full return address and daytime phone number; a note describing the problem you experienced; **a copy of your sales receipt or other proof of purchase to determine warranty status**. C.O.D. shipments cannot be accepted.

Return the appliance to: ATTN: Repair Department, Salton, Inc.
550 Business Center Drive, Mt. Prospect, Illinois 60056

For more information on Salton, Inc. products, visit our website:
<http://www.salton-maxim.com>, or email us at: salton@saltonusa.com

NOTE: Trademarks are used under license and should not be copied or used in any way without prior permission from the Licensor and/or Licensee

TM & © Dr. Seuss Enterprises, L.P. 1957. Licensed by Universal Studios Licensing, Inc.
All Rights Reserved.

Printed in China
©2000 Salton, Inc.

P/N 60582