


# RELAXOR<sup>®</sup>

*Perfect Touch<sup>™</sup> Air Massage System*

## ***Control Unit***

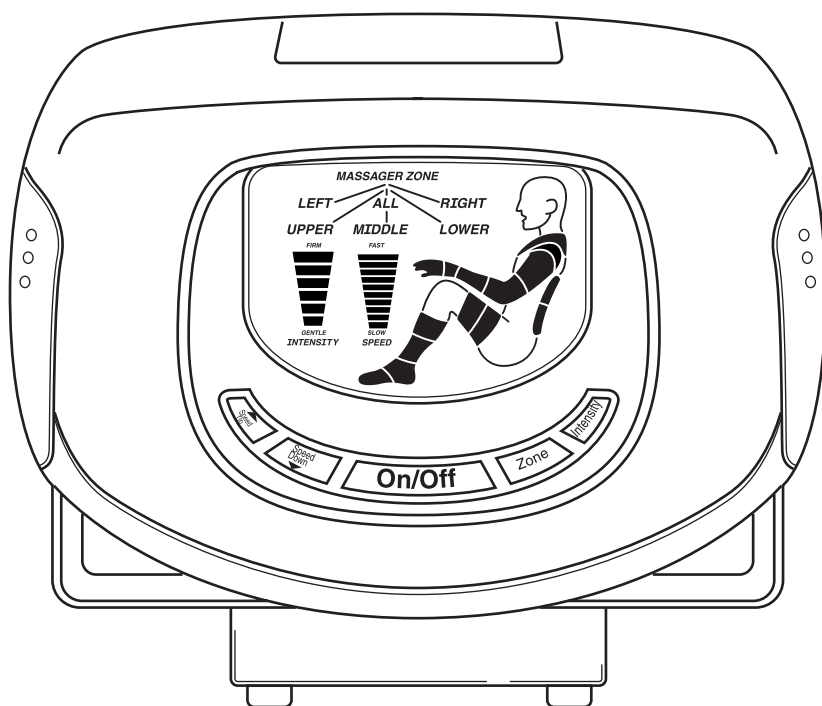

## **IMPORTANT SAFEGUARDS**

When using an electrical appliance, basic safety precautions should always be followed, including the following:

Read all instructions before using the appliance.

### **WARNING:**

This product is intended for use by people who enjoy good health

- This product is not recommended for people who use pace makers, have heart problems or vascular problems.
- This product is not recommended for people with spinal cord or back injuries.
- This product is not recommended for women who are pregnant.
- This product is not recommended if you are under physicians care, have a condition requiring the use of medical device, or any other condition that might affect your normal well being.
- Do not use this product over sensitive areas of the skin, or in the presence of poor circulation.
- Do not use this product on any part of the body that is swollen, infected, burned or inflamed, or where skin eruptions or sores are presented.
- Contact your physician prior to use.

**DANGER - This product contains electronic components**

**WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:**

1. Never leave this product unattended when plugged in. Always unplug this product from the electrical outlet immediately after using, before cleaning, and before putting on or taking off parts.
2. Use this product only for its intended use as described in this manual.
3. Never use pins or other metallic fasteners with this product.
4. Carefully examine the product before each use. Discontinue use of the product if it shows and sign of deterioration, such as checking, blistering, or cracking. Take product to the nearest authorized service center for examination and repair.
5. Keep Dry- Do not operate while bathing, in a shower, in or around water, or in a wet or moist condition. Do not use lotions, creams, gels, or similar skin products with this product

6. Keep this product out of the reach of children and pets.  
Close supervision is necessary when this product is used near children, invalids, or disabled persons. Use by children may be dangerous. Unattended use by incapacitated persons may be dangerous.
7. Use this product only for its intended use as described in the manual. Do not use attachments not recommended by the manufacturer.
8. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged in any manner, or dropped into water. Take product to the nearest authorized service center for examination and repair.
9. Do not reach for the product if it has fallen into water.  
Unplug immediately.
10. Do not carry this product by the supply cord or use the cord as a handle. Always use the built-in carrying handle.
11. Keep the cord away from heated surfaces.
12. Never block the air openings of the product or place it on a soft surface such as a bed, pillow, or couch, where air openings may be blocked. Keep the openings free of lint, hair, and the like.
13. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
14. Never use while sleeping.
15. Do not use outdoors.
16. Do not use where aerosol (spray) products are being used or oxygen is being administered.
17. To disconnect, turn all the controls to the OFF position, then remove the plug from the outlet.

**SAVE THESE INSTRUCTIONS**

**FOR HOUSEHOLD USE ONLY**

# CONGRATULATIONS!

You now own a top-quality Relaxor product – thoroughly tested to be certain it meets our own exacting quality control standards. Please take a few minutes to review the following important information before using this item.

## INDICATIONS FOR USE:


The *Perfect Touch*™ Air Massager System is indicated for temporary relief of minor muscle aches and pains and for temporary increase in circulation to the treated areas in people who are in good health. The *Perfect Touch*™ simulates kneading and stroking of tissues by using an inflatable garment.

## GETTING TO KNOW YOUR RELAXOR® *PERFECT TOUCH*™ MASSAGER

### HOW IT WORKS:

The *Perfect Touch*™ Air Massager stimulates kneading and stroking of tissue with the hands by using air pressure to alternately inflate and deflate a series of bladders or balloons within the Massage Garments. When the garments are worn on the body, they simulate the touch of the human hand by repeatedly pressing and releasing the body area—providing a soothing and comforting sensation, relief from everyday stress and tension, pampering and soothing. Further, your *Perfect Touch*™ Air Massager is designed as a massage SYSTEM. Many different Massage Garments are available. Each of the Massage Garments allows massaging of a different body area and is designed to work with your *Perfect Touch*™ control unit.

### Massage Sequence

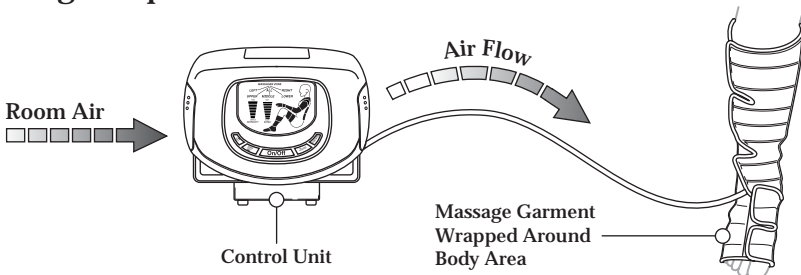


**RECOMMENDED MASSAGE TIME:**

Generally, massage should only last 10-15 minutes on any muscle of the body. For safety purposes, and to ensure comfort, the control unit will automatically shut off after 15 minutes of use. To restart the Air Massager, simply press the On/Off button on the front panel of the control unit. It is recommended that the Air Massager not be used for more than 30 minutes in one session. Excessive massaging can over stimulate muscles and nerves and result in an adverse effect. This product should only be used once or twice daily.

**FOR BEST RESULTS:**

Greatest sensation is achieved when worn directly against exposed skin. Can be worn over clothing. Clothing will lessen the effect if there is too much stimulation. If you feel that the Air Massager gives you too much stimulation, try placing a cloth or towel between you and the garment. This may act as a buffer to the massaging action. Wearing the Massage Garments over clothing is an effective way to buffer the massaging action. Massage should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your doctor.

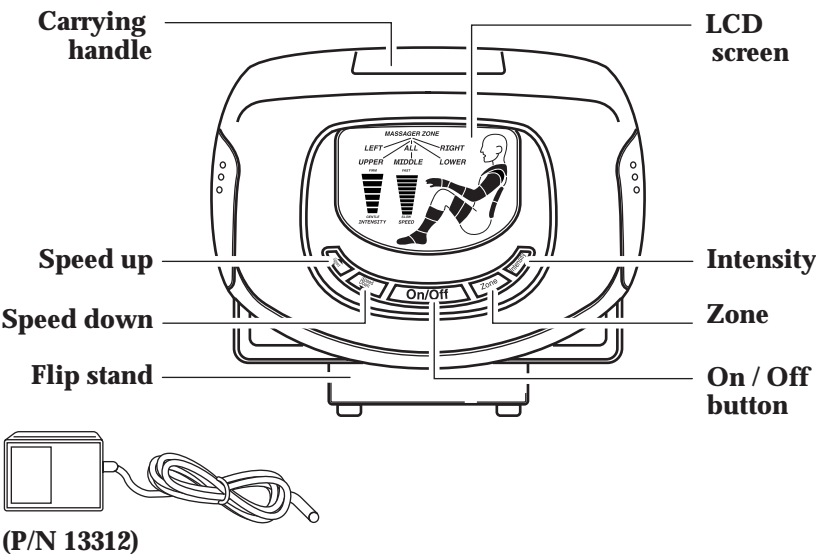


## INSTRUCTIONS FOR USE:

1. Plug the Adapter into the DC Jack socket on the back of the controller.
2. Plug the Adapter into a 120V AC electrical outlet

**Note:** Use only the adapter supplied with this product



3. Press the power On/Off switch on the back of the control unit to the “ON” position (this switch provides power to the control unit).
4. Next, Press the On/Off button on the control panel on the front of your unit to begin your massage.
5. The red LED light will illuminate to indicate that the massager is operating.

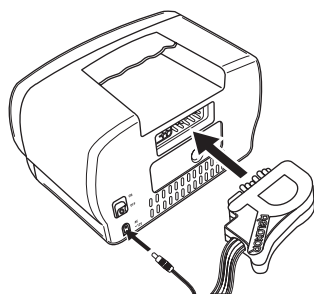

### Attaching the Massage Garments:

To connect your Massage Garment, follow the below instructions

1. Locate the “R” shaped plug at the end of the tubing attached to your Massage Garment.
2. Locate the plug receptor on the back of the unit.
3. Line up the plug and gently press plug into the receptor.
4. The *Perfect Touch™* Air Massager will automatically recognize that you have attached the Massage Garment you have attached.
5. The portion of the human figure on the LCD screen will begin to flash if garment is plugged in correctly.

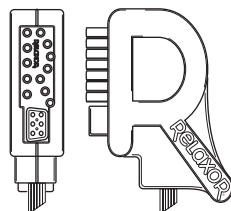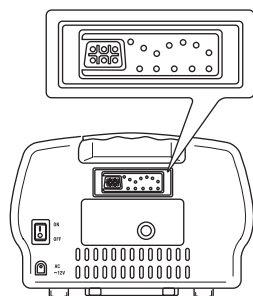

## TO WEAR GARMENT:

Control Unit and two massage garments are included with the Control Unit in this set. The Leg and Foot garment, and the Thigh garment will help you to begin building your Perfect Touch Air Massage System. Please refer to individual garment manuals for more specific instructions on these

**Note:** The Massage Garment will fit either the left or right body part being massaged (except for the Hand Massage Garment and the Feet Massage Garment)


1. Undo all the Velcro tabs on the garment and open completely flat. Hold the garment by its outermost corners of the top edge and position body part in the Massage Garment.
2. Refer to your Massage Garment Manual for more specific instructions on wearing that particular Massage Garment.

**Note:** ONE SIZE FITS MOST. The garment should be snug at all points—not too tight, not too loose—readjust individual tabs as necessary. Properly worn the garment should look as shown.

**TIP:** The better the garment fits, the better the massaging action will be. For best results, close the Velcro tabs in the order shown— from the toes to the knee.

**TIP:** Although the garment can be worn over clothing, for the maximum effect wear directly against the skin.

**TIP:** Try to engage as much of the Velcro tabs as possible - at least 3/4" if possible. This will ensure that the garment will stay closed at the higher intensity levels.

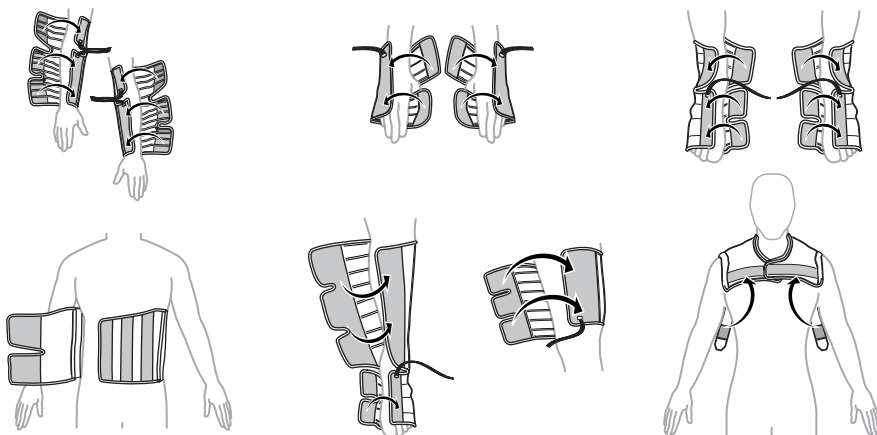


## OPERATION:

1. Once Massage Garment is secured around area to be massaged, press the On/Off button on the front control panel.
2. Garment will begin to inflate. Allow a couple of inflation cycles through the garment before selecting Speed, Intensity, or Zone.

### Functions:

#### SPEED:

The speed setting determines how fast or slow the inflation cycle will take. A higher speed will cause each bladder to fill up a quicker rate, while lower speed will cause each bladder to fill at a slower rate. There are 10 SPEED settings from SLOW to FAST.

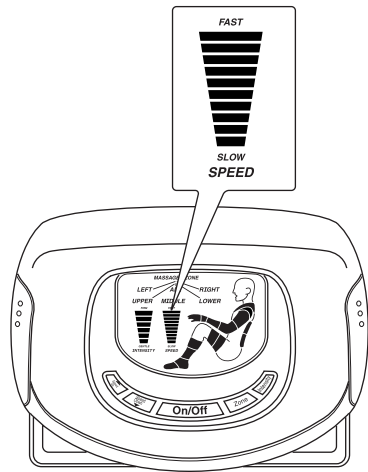

### To adjust speed:


1. To increase speed, press the Speed Up button. On the LCD screen, you will see more bars appear in the speed gauge to indicate a higher speed.
2. To decrease speed, press the Speed Down button. On the LCD screen, you will see bars disappear on the speed gauge to indicate a lower speed



## OPERATION: continued

### INTENSITY:

The intensity setting determines how much air pressure (psi's) the control unit pushes into the Massage Garment. A higher intensity will cause each bladder to fill up with more air pressure, while lower speed will cause each bladder to fill with less air pressure. There are 6 levels of INTENSITY from GENTLE to FIRM.

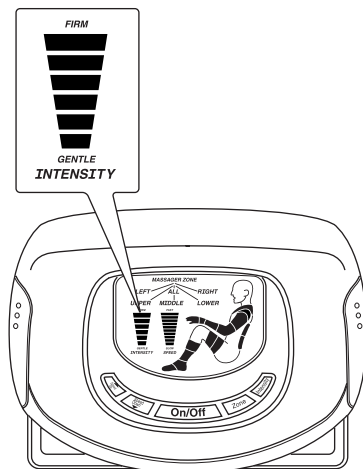

### To adjust intensity:

1. To increase speed, press the Intensity button. On the LCD screen, you will see more bars appear in the intensity gauge to indicate a higher intensity.
2. To decrease speed, press the Intensity button. On the LCD screen, you will see bars disappear on the intensity gauge to indicate a lower intensity.

**Note:** at higher intensities there may be an audible loosening of the Velcro. This is usual. If any of the Velcro tabs open at any time, simply stop the massager and adjust the Leg and foot garment for a better fit and closure, then restart.


## ZONE:

The ZONE feature allows you to concentrate the massage on a particular area. Each Massage Garment has specific zones available. Refer to your Massage Garment Manual for zones available for that particular body part.

### To adjust the ZONE:

1. Press ZONE button repeatedly until you have chosen the exact body location you want continuously massaged.
2. The corresponding massage zones will be displayed to indicate the area you have targeted.
3. The LCD display will flash in the area you have targeted, and the words will also indicate the ZONE you have selected.
4. The massaging action will remain in the areas you have selected until you press the ZONE button again, which will change the massage to the next Massage Zone.
5. If you want to massage the entire body area, repeatedly press the ZONE button until ALL is selected.

**Note:**The *Perfect Touch*<sup>TM</sup> Air Massager defaults to a power-on setting approximating the mid-point of the SPEED and INTENSITY settings and automatically zones to ALL.

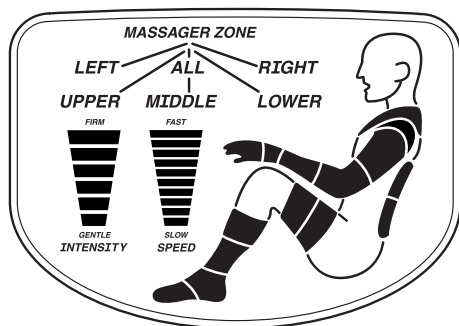


## TO REMOVE THE GARMENT:

To remove garment:


1. Stop the massage session by pressing the On/Off button on the front control panel. Then simply pull the Velcro tabs to release the garment. Please fold the garment carefully and coil the airline carefully for storage.

### Tilt Stand:

*Tilt stand shown in closed position*

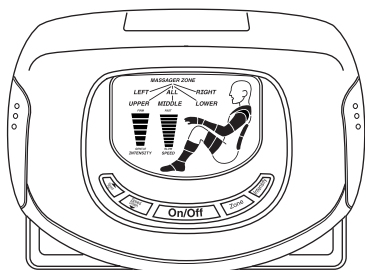


*Tilt stand shown in open position*

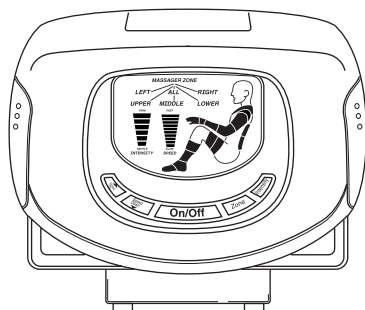

The *Perfect Touch*™ Air Massager has an integral tilt stand feature to allow for ease of view for either table-top or floor placement. Simply pull the stand down to use; push back up to store. Use the stand as necessary for optimal viewing and positioning of the LCD screen and touch controls.

### Carrying handle:

The *Perfect Touch*™ Air Massager has a built-in carrying handle. Always carry the unit by this handle. Never use the power cord or Massage Garments' cord to carry the unit; damage will occur.

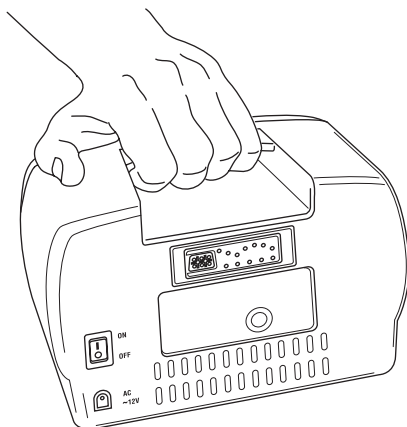

## MASSAGE GARMENTS:


A variety of Massage Garments are available for massaging different body parts. All Massage Garments are 100% fully compatible with your *Perfect Touch™* Air Massager control unit - simply plug in and turn on!

**Note: ONE SIZE FITS MOST.**


**AM8AR**



**AM8FEET**

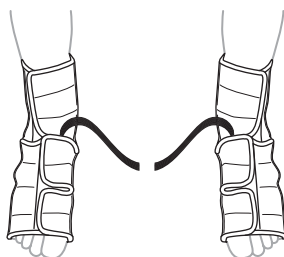

**AM8HAND**

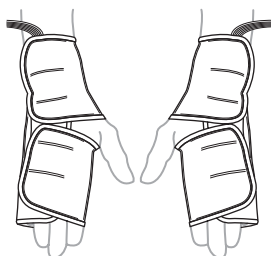


**AM8LF**




**AM8THIGH**



**AM8LB**

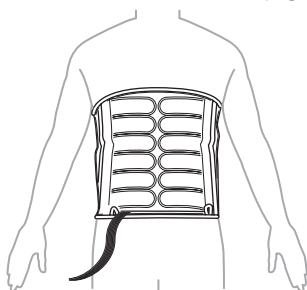

**AM8NS**

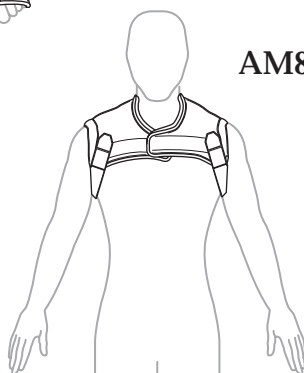

**TO CLEAN: NEVER IMMERSE GARMENT IN WATER OR ANY OTHER LIQUID.** Use a soft, slightly damp cloth or sponge to wipe the garment surface. **DO NOT SATURATE.** Allow adequate time to dry before using. Never use abrasive cleaners or solvents, as they may damage the surface.

### **POLARIZED PLUG**

This appliance has a polarized plug (one blade wider than the other). As a safety feature, this plug fits into a polarized outlet one way. If the plug does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

**TO STORE:** Unplug garment before storing. The garment can be folded for space saving storage. Coil the air line. Avoid contact with sharp edges or pointed objects that might cut or puncture the fabric surface. Use care not to bend the airline too severely to prevent kinking. Store at room temperature.

Notes:

---

---

---

---

---

---

---

---

---

---

---

---

---

---

---

# LIMITED ONE YEAR WARRANTY

**Warranty:** This Relaxor® product is warranted by Salton, Inc. to be free from defects in materials or workmanship for a period of (1) year from the original purchase date. This product warranty covers only the original consumer purchaser of the product.

**Warranty Coverage:** This warranty is void if the product has been damaged by accident in shipment, unreasonable use, misuse, neglect, improper service, commercial use, repairs by unauthorized personnel, normal wear and tear, improper assembly, installation or maintenance abuse or other causes not arising out of defects in materials or workmanship. This warranty is effective only if the product is purchased and operated in the USA, and does not extend to any units which have been used in violation of written instructions furnished with the product or to units which have been altered or modified or to damage to products or parts thereof which have had the serial number removed, altered, defaced or rendered illegible.

**Implied Warranties:** ANY IMPLIED WARRANTIES WHICH THE PURCHASER MAY HAVE ARE LIMITED IN DURATION TO ONE (1) YEAR FROM THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

**Warranty Performance:** During the above one-year warranty period, a product with a defect will be either repaired or replaced with a reconditioned comparable model (at our option) when the product is returned to the Repair Center, or the purchase price refunded. The repaired or replacement product will be in warranty for the balance of the one-year warranty period and an additional one-month period. No charge will be made for such repair or replacement.

**Service and Repair:** Should the appliance malfunction, you should first call toll-free 1 (800) 934-7455 between the hours of 9:00 am and 5:00 pm Central Standard Time and ask for CONSUMER SERVICE stating that you are a consumer with a problem. Please refer to model number 61131 when you call.

**In-Warranty Service (USA):** For an appliance covered under the warranty period, no charge is made for service or postage. Call for return authorization (1 800 934-7455).

**Out-of-Warranty Service:** A flat rate charge by model is made for out-of-warranty service. Include \$20.00 (U.S.) for return shipping and handling. We will notify you by mail of the amount of the charge for service and require you to pay in advance for the repair or replacement.

**For Products Purchased in the USA, but Used in Canada:** You may return the product insured, packaged with sufficient protection, and postage and insurance prepaid to the USA address listed below. Please note that all customs duty / brokerage fees, if any, must be paid by you and we will require you to pay the cost of customs duty / brokerage fees to us in advance of our performing any service.

**Risk During Shipment:** We cannot assume responsibility for loss or damage during incoming shipment. For your protection, carefully package the product for shipment and insure it with the carrier. Be sure to enclose the following items with your appliance: any accessories related to your problem, your full return address and daytime phone number, a note describing the problem you experienced, a copy of your sales receipt or other proof of purchase to determine warranty status. C.O.D. shipments cannot be accepted.

**To return the appliance, ship to:**

Attn: Repair Center  
Salton, Inc.  
708 South Missouri Street  
Macon, MO 63552

**To contact us, please write to, call, or email:**

Salton, Inc.  
P.O. Box 1637  
Columbia, MO 65205-1637  
1(800) 934-7455  
Email: Salton@Saltonusa.com

**Limitation of Remedies:** No representative or person is authorized to assume for Salton, Inc. any other liability in connection with the sale of our products. There shall be no claims for defects or failure of performance or product failure under any theory of tort, contract or commercial law including, but not limited to negligence, gross negligence, strict liability, breach of warranty and breach of contract. Repair, replacement or refund shall be the sole remedy of the purchaser under this warranty, and in no event shall Salton be liable for any incidental or consequential damages, losses or expenses. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

**Legal Rights:** This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

**For more information on Salton products:** Visit our website: [www.eSalton.com](http://www.eSalton.com)