


AM8FEET

RELAXOR[®]

Perfect Touch™ Air Massage System

Feet Garment

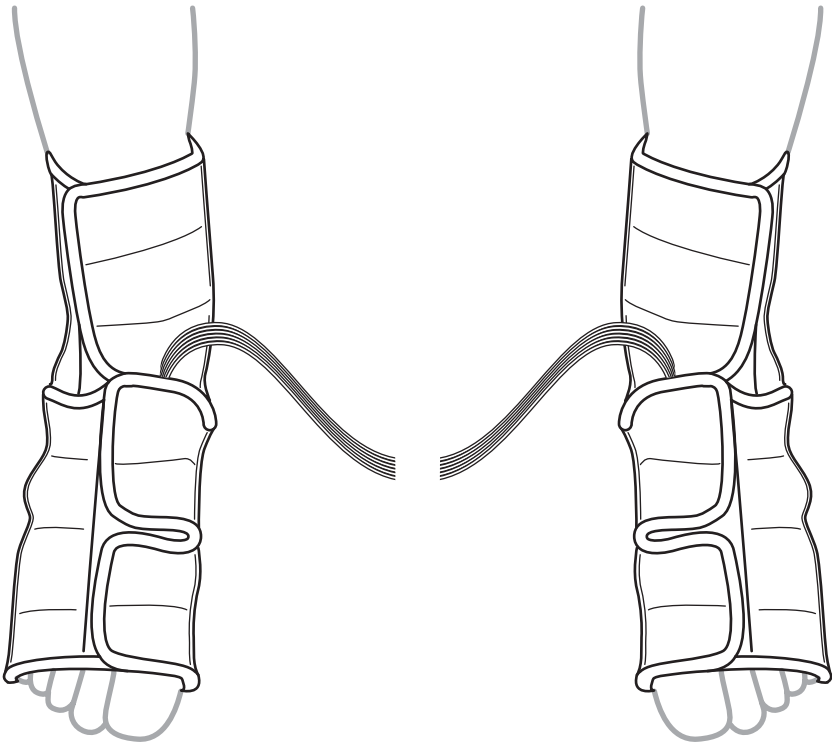

IMPORTANT SAFEGUARDS

Read all instructions before using the appliance. When using an electrical appliance, basic safety precautions should always be followed, including the following:

DANGER: To reduce the risk of electrical shock:

- Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- Do not reach for massage garment that has fallen into water. Unplug it immediately.
- Do not use while bathing or showering.
- Do not place or store massage garment where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or any other liquid.
- Do not attempt to plug or unplug massage garment while you are standing in water.
- Do not use over wet skin or wet clothing.

WARNING:

WARNING:

This product is intended for use by people who enjoy good health

- This product is not recommended for people who use pace makers, have heart problems or vascular problems.
- This product is not recommended for people with spinal cord or back injuries.
- This product is not recommended for women who are pregnant.
- This product is not recommended if you are under physicians care, have a condition requiring the use of medical device, or any other condition that might affect your normal well being.
- Do not use this product over sensitive areas of the skin, or in the presence of poor circulation.
- Do not use this product on any part of the body that is swollen, infected, burned or inflamed, or where skin eruptions or sores are presented.
- Contact your physician prior to use.

! WARNING: To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Massage garment must never be immersed in water.
2. Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
3. Use this product only for its intended use as described in the manual. Use only those attachments recommended by the manufacturer.
4. Never use massage garment while sleeping or fall asleep while using the massager.
5. Ensure hands are dry when connecting the massage garment to the control unit.
6. Do not place or drop into water, or any other liquid.
7. Never use the massage garment if it has a damaged cord or plug, if it is not working properly, or if has been damaged or dropped in water. Please contact consumer service for return authorization.
8. Keep cord away from heated surfaces.
9. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
10. This massage garment is designed for household use only.

CAUTION: Please read all instructions carefully before operating.

- If you have concerns regarding your health, consult your doctor before using the massage garment.
- Individuals with pacemakers should consult a physician before using this product
- Do not use this product for more than 30 continuous minutes.
- Never use massage garment on swollen or inflamed areas or on skin eruptions.
- Massage garment should not be used by children without close adult supervision.
- In case of pregnancy, diabetes, or illness, consult your doctor before using the massage garment
- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult a doctor.
- Massage should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.

SAVE THESE INSTRUCTIONS

FOR HOUSEHOLD USE ONLY

CONGRATULATIONS!


You now own a top-quality Relaxor® product- thoroughly tested to be certain it meets our own exacting quality control standards. Please take a few minutes to review the following important information before using this item.

RECOMMENDED MESSAGE TIME:

Generally, massage should only last 10-15 minutes on any part of the body. For safety purposes, and to ensure comfort, the control unit will automatically shut off after 15 minutes of use. To restart the air massager, simply press the On/Off button on the front control panel. It is recommended that the Air Massager not be used for more than 30 minutes in one session. Excessive massaging can over-stimulate muscles and nerves and result in an adverse effect. This product should only be used once or twice daily.

FOR BEST RESULTS:

Greatest sensation is achieved when worn directly against exposed, dry skin. However, it can be worn over clothing. Clothing will lessen the effect if there is too much stimulation. If you feel that the Air Massager gives you too much stimulation, try placing a cloth or towel between you and the garment. This may act as a buffer to the massaging action. Wearing the Massage Garments over clothing is an effective way to lessen the massaging action. Massage should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your doctor.

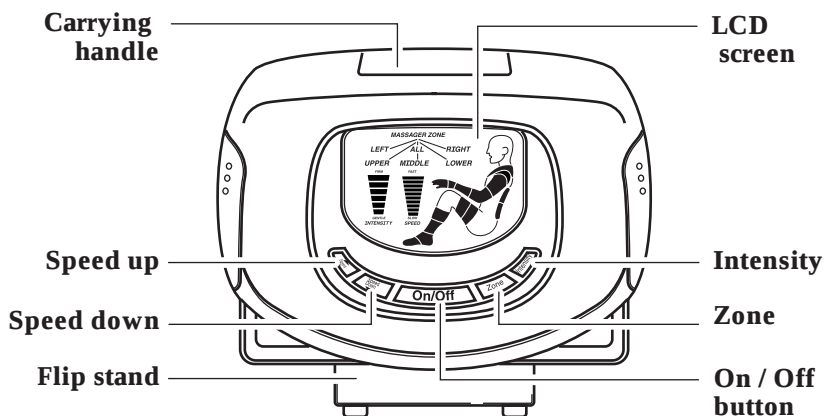


INSTRUCTIONS FOR USE:

1. Plug the Adapter into the DC Jack socket on the back of the controller.
2. Plug the Adapter into a 120V AC electrical outlet

Note: Use only the adapter supplied with this product


3. Press the power button on the back of the unit to the on position.
4. Press the On/Off button on the front panel of your control unit to begin your massage.
5. The red LED light will illuminate to indicate that the massager is operating.

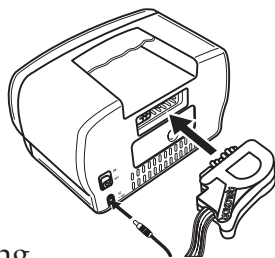

ATTACHING THE FOOT MESSAGE GARMENT


To connect your Foot Massage garment, follow the below instructions

1. Locate the “R” shaped plug at the end of the tubing attached to your Foot massage garment.
2. Locate the plug receptor on the back of the unit.
3. Line up the plug and gently press plug into the receptor. (The word “Relaxor®” should be face up)
4. The Perfect Touch™ Air Massager will automatically recognize that you have attached the Foot Garment.
5. The Foot portion of the human figure on the LCD screen will begin to flash if garment is plugged in correctly.

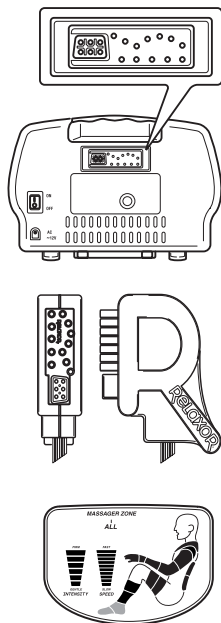

TO WEAR GARMENT:

1. Undo all the Velcro tabs on the garment and open completely flat. Hold the garment by its outermost corners of the top edge and position Feet in the garment.
2. Place your foot on the lower "foot" section of the garment with your heel snugly inserted in the round cut-out opening. Wrap the straight edge (not the tabbed edge) over the foot, hold in position and then close the two lower tabs around the foot - aim for a snug fit..
3. Repeat for the ankle area - wrap and hold the straight edge while closing the single tab around the ankle

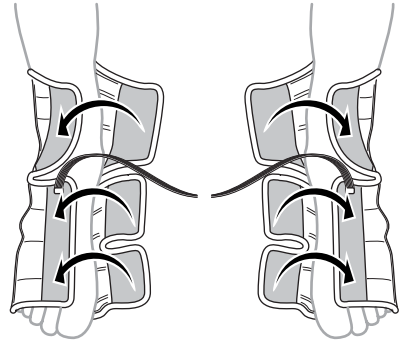


Note: ONE SIZE FITS MOST. The garment should be snug at all points - not too tight, not too loose. Readjust individual tabs as necessary. Properly worn the garment should look as shown.

TIP: The better the garment fits, the better the massaging action will be.

TIP: Although the garment can be worn over clothing, for the maximum effect wear directly against the skin.

TIP: Try to engage as much of the Velcro tabs as possible - at least 3/4" if possible. This will ensure that the garment will stay closed at the higher intensity levels.

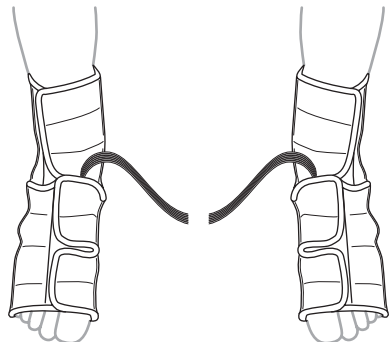

OPERATION:

1. Once massage garment is secured around area to be massaged, press the On/Off button on the front control panel.
2. Garment will begin to inflate. Allow a couple of inflation cycles through the garment before selecting Speed, Intensity, or Zone.

Perfect Touch Speed / Intensity settings -- Feet Garment											
Speed Settings	1	2	3	4	5	6	7	8	9	10	
Inflation (seconds)	6.0	5.5	4.87	4.25	3.62	3.0	2.37	1.75	1.12	0.50	
Intensity Settings	1		2		3		4		5		6
Garment Air Pressure	80 (mmHg)		104 (mmHg)		128 (mmHg)		152 (mmHg)		176 (mmHg)		200 (mmHg)


Functions:

SPEED:

The speed setting determines how fast or slow the inflation cycle will take. A higher speed will cause each bladder to fill up at a quicker rate, while a lower speed will cause each bladder to fill at a slower rate. There are 10 SPEED settings from SLOW to FAST.

To adjust speed:

1. To increase speed, press the Speed Up button. On the LCD screen, you will see more bars appear in the speed gauge to indicate a higher speed.
2. To decrease speed, press the Speed Down button. On the LCD screen, you will see bars disappear on the speed gauge to indicate a lower speed.

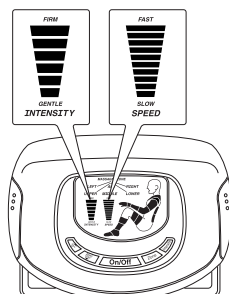

INTENSITY:

The intensity setting determines how much air pressure the control unit pushes into the massage garment. A higher intensity will cause each bladder to fill up with more air pressure, while lower intensity will cause each bladder to fill with less air pressure. There are 6 levels of INTENSITY from GENTLE to FIRM.

To adjust intensity:

1. To increase the massaging action, press the Intensity button. On the LCD screen, you will see more bars appear in the intensity gauge to indicate a higher intensity.
2. To decrease massaging action, press the Intensity button. On the LCD screen, you will see bars disappear on the intensity gauge to indicate a lower intensity.

Note: at higher intensities, and/or lower speeds, there may be an audible loosening of the Velcro. This is usual. If any of the Velcro tabs open at any time, simply stop the massager and adjust the Leg and Foot garment for a better fit and closure, then restart.


ZONE:

The ZONE feature allows you to concentrate the massage on a particular area. The Foot attachment has 5 available zones.


To adjust the ZONE:

1. Press ZONE button repeatedly until you have chosen the exact body location you want continuously massaged.
2. The corresponding massage zones will be displayed to indicate the area you have targeted.
3. The LCD display will flash in the area you have targeted, and the words will also indicate the ZONE you have selected.
4. The massaging action will remain in the areas you have selected until you press the ZONE button again, which will change the massage to the next Massage Zone.
5. If you want to massage the entire body part, repeatedly press the ZONE button until ALL is selected.
6. See the figure below for Zones available for your Foot Garment.


ALL
the entire garment

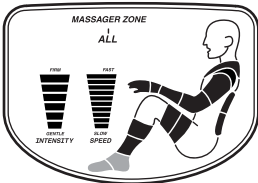


**LEFT + RIGHT
UPPER**
left & right upper foot area

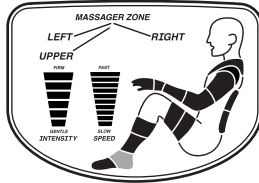


**LEFT + RIGHT
LOWER**
left & right lower foot area

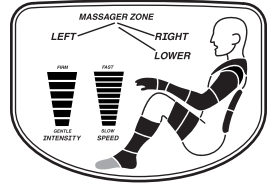

**LEFT
UPPER + LOWER**
left upper & lower foot area

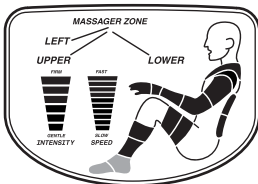

**RIGHT
UPPER + LOWER**
right upper & lower foot area

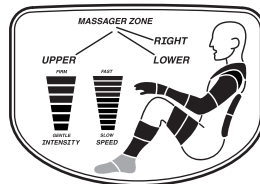

Note: The Perfect Touch™ Air Massager defaults to a power-on setting approximating the mid-point of the SPEED and INTENSITY settings and automatically zones to ALL.

To remove the garment:

1. Stop the massage session by pressing the On/Off button on the front control panel. Then simply pull the Velcro tabs to release the garment. Please fold the garment carefully and coil the airline carefully for storage.

TO CLEAN:

NEVER IMMERSE GARMENT IN WATER OR ANY OTHER LIQUID. Use a soft, slightly damp cloth or sponge to wipe the garment surface. **DO NOT SATURATE.** Allow adequate time to dry before using. Never use abrasive cleaners or solvents, as they may damage the materials.

TO STORE:


Unplug garment before storing. The garment can be folded for space saving storage. Coil the air line. Avoid contact with sharp edges or pointed objects that might cut or puncture the fabric surface. Use care not to bend the airline too severely to prevent kinking. Store indoors, at room temperature.

INSTRUCTIONS FOR USE:

1. Plug the Adapter into the DC Jack socket on the back of the controller.
2. Plug the Adapter into a 120V AC electrical outlet

Note: Use only the adapter supplied with this product


3. Press the power button on the back of the unit to the on position.
4. Press the On/Off button on the front panel of your control unit to begin your massage.
5. The red LED light will illuminate to indicate that the massager is operating.

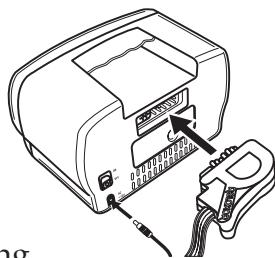

ATTACHING THE FOOT MESSAGE GARMENT

To connect your Foot Massage garment, follow the below instructions

1. Locate the “R” shaped plug at the end of the tubing attached to your Foot massage garment.
2. Locate the plug receptor on the back of the unit.
3. Line up the plug and gently press plug into the receptor. (The word “Relaxor®” should be face up)
4. The Perfect Touch™ Air Massager will automatically recognize that you have attached the Foot Garment.
5. The Foot portion of the human figure on the LCD screen will begin to flash if garment is plugged in correctly.

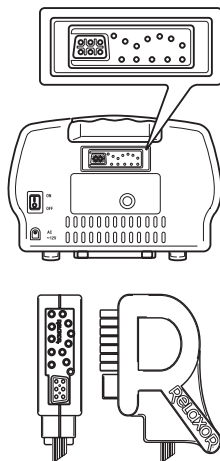

ADDITIONAL MASSAGE GARMENTS

Create the ultimate massage system
with these additional massage garments

Note: ONE SIZE FITS MOST.

AM8LF




AM8AR




AM8THIGH




AM8HAND

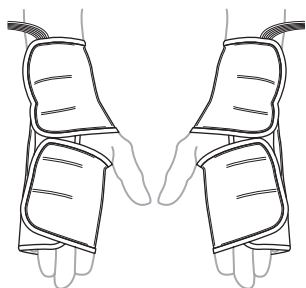

AM8NS

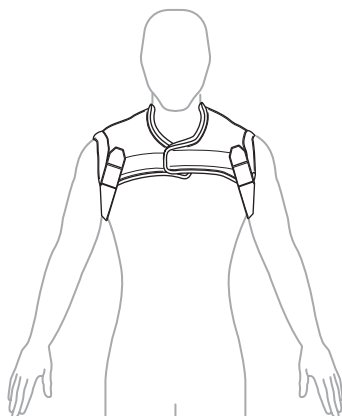

AM8LB

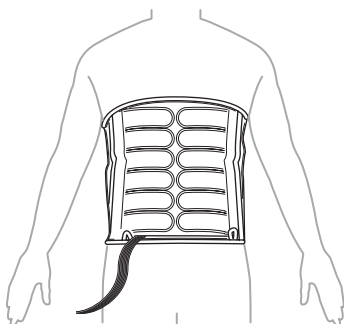

Notes:

LIMITED ONE YEAR WARRANTY

Warranty: This Relaxor® product is warranted by Salton, Inc. to be free from defects in materials or workmanship for a period of (1) year from the original purchase date. This product warranty covers only the original consumer purchaser of the product.

Warranty Coverage: This warranty is void if the product has been damaged by accident in shipment, unreasonable use, misuse, neglect, improper service, commercial use, repairs by unauthorized personnel, normal wear and tear, improper assembly, installation or maintenance abuse or other causes not arising out of defects in materials or workmanship. This warranty is effective only if the product is purchased and operated in the USA, and does not extend to any units which have been used in violation of written instructions furnished with the product or to units which have been altered or modified or to damage to products or parts thereof which have had the serial number removed, altered, defaced or rendered illegible.

Implied Warranties: ANY IMPLIED WARRANTIES WHICH THE PURCHASER MAY HAVE ARE LIMITED IN DURATION TO ONE (1) YEAR FROM THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

Warranty Performance: During the above one-year warranty period, a product with a defect will be either repaired or replaced with a reconditioned comparable model (at our option) when the product is returned to the Repair Center, or the purchase price refunded. The repaired or replacement product will be in warranty for the balance of the one-year warranty period and an additional one-month period. No charge will be made for such repair or replacement.

Service and Repair: Should the appliance malfunction, you should first call toll-free 1 (800) 934-7455 between the hours of 9:00 am and 5:00 pm Central Standard Time and ask for CONSUMER SERVICE stating that you are a consumer with a problem. Please refer to model number 61126 when you call.

In-Warranty Service (USA): For an appliance covered under the warranty period, no charge is made for service or postage. Call for return authorization (1 800 934-7455).

Out-of-Warranty Service: A flat rate charge by model is made for out-of-warranty service. Include \$20.00 (U.S.) for return shipping and handling. We will notify you by mail of the amount of the charge for service and require you to pay in advance for the repair or replacement.

For Products Purchased in the USA, but Used in Canada: You may return the product insured, packaged with sufficient protection, and postage and insurance prepaid to the USA address listed below. Please note that all customs duty / brokerage fees, if any, must be paid by you and we will require you to pay the cost of customs duty / brokerage fees to us in advance of our performing any service.

Risk During Shipment: We cannot assume responsibility for loss or damage during incoming shipment. For your protection, carefully package the product for shipment and insure it with the carrier. Be sure to enclose the following items with your appliance: any accessories related to your problem, your full return address and daytime phone number, a note describing the problem you experienced, a copy of your sales receipt or other proof of purchase to determine warranty status. C.O.D. shipments cannot be accepted.

To return the appliance, ship to:

Attn: Repair Center
Salton, Inc.
708 South Missouri Street
Macon, MO 63552

To contact us, please write to, call, or email:

Salton, Inc.
P.O. Box 1637
Columbia, MO 65205-1637
1(800) 934-7455
Email: Salton@Saltonusa.com

Limitation of Remedies: No representative or person is authorized to assume for Salton, Inc. any other liability in connection with the sale of our products. There shall be no claims for defects or failure of performance or product failure under any theory of tort, contract or commercial law including, but not limited to negligence, gross negligence, strict liability, breach of warranty and breach of contract. Repair, replacement or refund shall be the sole remedy of the purchaser under this warranty, and in no event shall Salton be liable for any incidental or consequential damages, losses or expenses. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

Legal Rights: This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

For more information on Salton products: Visit our website: www.eSalton.com

© 2002 Salton, Inc.
All Rights Reserved

P/N 61126
Printed in China