


RELAXOR®

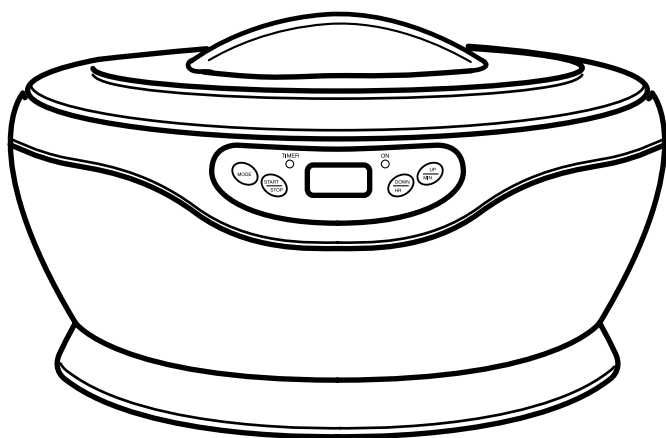

Electronic Paraffin Wax Bath

OWNER'S MANUAL &
OPERATING INSTRUCTIONS

Model PWB8FE

SPA THERAPY™

IMPORTANT SAFEGUARDS

When using electrical appliances, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGER: *To reduce the risk of electrocution.*

- 1. Always unplug this appliance immediately after using.*
- 2. Do not use while bathing.*
- 3. Do not place or store this appliance where it can fall or be pulled into a tub or sink.*
- 4. Do not place or drop into water or any other liquid.*
- 5. Do not reach for an appliance that has fallen into water. Unplug it immediately.*

WARNING: *To reduce the risk of burns, electrocution, fire or injury to persons:*

- 1. An appliance should not be left unattended when plugged in.*
- 2. Close supervision is necessary when this unit is used by, on or near children or invalids.*
- 3. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by Relaxor(.*
- 4. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged in any manner or if it has been dropped into water. Refer to the warranty section of this manual.*
- 5. Keep the Power Cord away from heated surfaces.*
- 6. Do not carry the appliance by the cord, or use the cord as a handle.*
- 7. Do not operate the unit without the paraffin wax.*
- 8. Never place the unit on a soft surface such as a bed or couch. Always place unit on a stable, level surface.*
- 9. Never use while sleepy or drowsy.*
- 10. Unplug the unit before filling. Do not overfill.*

SAVE THESE INSTRUCTIONS FOR HOUSEHOLD USE ONLY

ADDITIONAL IMPORTANT SAFEGUARDS

1. Do not use an extension cord with this appliance.
2. Use only recommended paraffin wax in this appliance.
3. Do not use paraffin wax near eyes or other heat-sensitive areas, or on any skin injury, eruption or surface growth.
4. Never use sharp instruments to remove solidified paraffin wax from the unit.
5. Test hot paraffin on a small area of skin for temperature and possible skin reactions before using it. If appliance is to be used on persons susceptible to allergic reactions or suffering from chronic circulatory problems it is recommended that they consult a physician before using unit.
6. Do not use paraffin wax on open cuts, wounds, warts, inflamed skin or where sensation in extremities is decreased.
7. Discontinue use and consult your physician if pain occurs from continued paraffin treatments.
8. Do not put water into the unit.
9. Do not use while in bed.
10. Use this appliance only on a level, heat-resistant surface.
11. Unplug unit when not in use. Allow unit to cool before handling, cleaning or storing the appliance.
12. **Do not put more than five pounds of wax in the unit at any one time.**

POLARIZED PLUG

This appliance has a polarized plug (one blade is wider than the other). This plug is intended to fit into a polarized outlet only one way. If the plug does not fit fully in the outlet reverse the plug. If it still does not fit contact a qualified electrician. Do not attempt to modify the plug in any way.


PLASTICIZER WARNING

CAUTION: To prevent Plasticizers from migrating from the finish of the counter top or table top or other furniture, place **NON-PLASTIC** coasters or place mats between the appliance and the finish of the counter top or table top.

Failure to do so may cause the finish to darken, permanent blemishes may occur or stains can appear.

SETTING UP YOUR WAX BATH

1. *Unpack unit and all accessories. Wipe all surfaces with a paper towel or a clean cloth.*
2. *Unwrap cord and run it through the Cord Slot.*
3. *Place unit on a sturdy, level surface.*
4. *Remove Liner Lid and put Liners in the Liner Storage Area; replace Liner Lid using the hook to secure.*
5. *Remove Lid and pour four lbs. of paraffin wax into the Wax Chamber; replace Lid.*

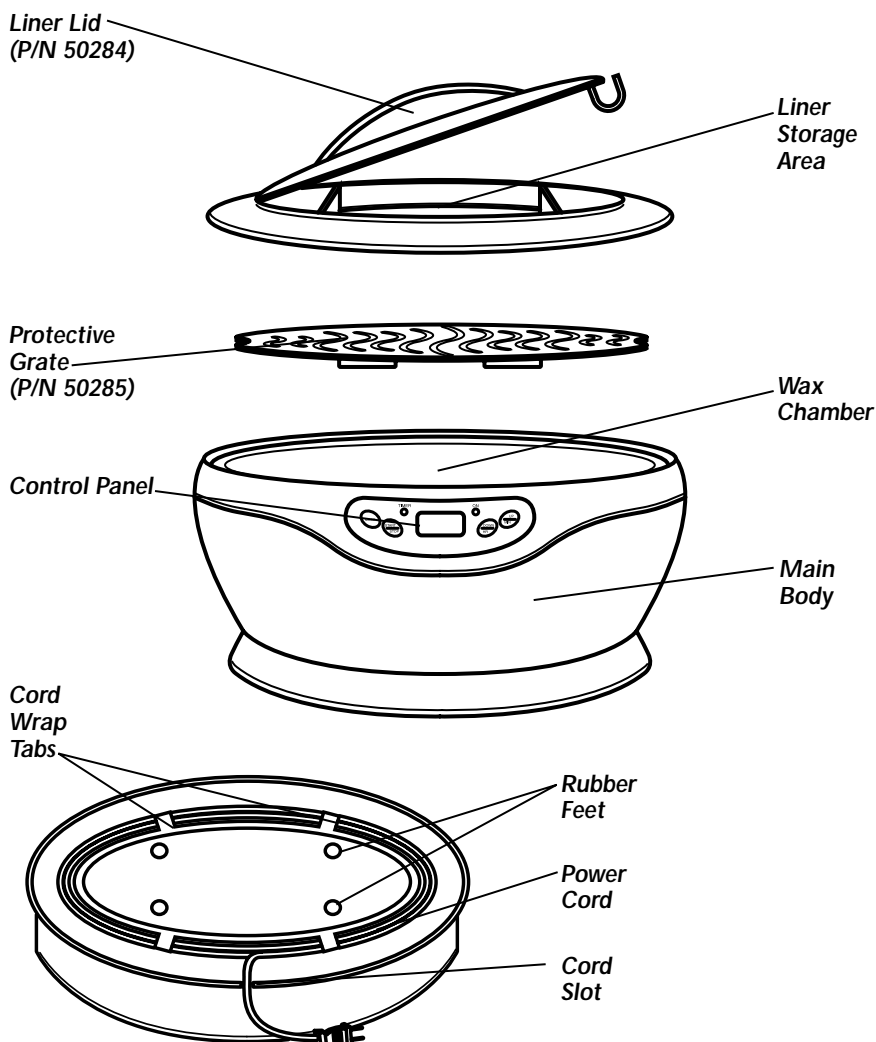

SETTING THE CLOCK

The Clock and Timer features allow you to program the unit to begin the melting process automatically. If you do not wish to use these features, and want to start the melting process immediately, skip to the SETTING THE TEMPERATURE section.

1. Plug unit into standard household 120V outlet.
2. 12:00 will flash on the Display.
3. Press the DOWN/HR Button to advance the hour.
4. Press the UP/MIN Button to advance the minute.
5. Once time is correct, press the START/STOP Button to set.



SETTING THE TIMER

1. Once the time has been set, press the MODE Button twice to bring up the Timer Display. Before setting, the Timer Display will read 0:00.
2. Press the DOWN/HR Button to adjust set the Timer hour.
3. Press the UP/MIN Button to adjust set the Timer minute.
4. Once you have selected the time you would like melting to begin, press the START/STOP Button.
5. The Timer is now set; the Timer LED will illuminate. The unit will automatically begin melting the wax at the selected time.

NOTE: The Clock and Timer will only remain set as long as the unit is plugged in. They will have to be reset each time you wish to unplug the unit.

SETTING THE TEMPERATURE

1. Press the Mode Button until the Display shows the temperature setting. The Display will show the default temperature of 131°F.
2. Press the UP/MIN and DOWN/HR Buttons until the unit displays between 122°F - 133°F.

NOTE: Wax can take a long time to change temperatures once it has melted. It is recommended that you select the desired temperature before melting the wax.

NOTE: The wax will take longer to melt at lower temperatures.

INSTRUCTIONS FOR USE

This unit is designed to allow for dipping hands, elbows and feet. Do not dip any other body parts.

1. If the Timer Function is not being used, press the START/STOP Button to activate unit. If the Clock has not been set, you will have to press the START/STOP Button twice to activate unit. The green ON Light will flash to indicate that wax is melting.
2. After approximately 4 hours, the green Light will stop flashing and remain steady. This indicates that the unit is ready for use.
3. Place the Protective Grate in unit; allow it to sink to the floor of the Wax Chamber.
4. Remove all jewelry from areas you intend to dip. Examine skin for cuts, infections or rashes. If any are present do not use the Paraffin Wax Bath.

5. Wash hands thoroughly. If desired, apply lotion to areas you intend to dip. This will allow for easier removal of the paraffin wax.

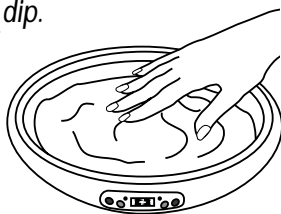



Illustration A

6. Dip hand, foot or elbow in the wax. **WAX WILL BE HOT! DO NOT LEAVE ANY BODY PART IN THE WAX FOR MORE THAN A FEW SECONDS.** Repeat dipping five to nine times to ensure good coverage. For best coverage, spread fingers or toes slightly to allow the wax to coat all areas evenly. See Illustration A.



Illustration B
Plastic Liners
(P/N 50252)

7. Slip hand or foot into Plastic Liner. See Illustration B.
8. Repeat steps 5-8 with other hand or foot.

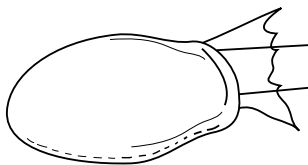

9. Slip hands or feet into Mitts or Booties. See Illustration C.

NOTE: Do not walk with Plastic Liners and/or Booties on your feet, this may cause you to slip and fall.

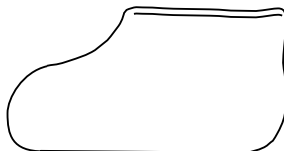


Illustration C
Mitts (P/N 50286)
Booties (P/N 50287)

HELPFUL TIPS

1. If the surface of the paraffin hardens (forms a skin), raise the temperature by pressing the UP/MIN Button while in the Temperature Mode.
2. Use only paraffin wax recommended by Relaxor® in the Paraffin Wax Bath. Do not attempt to use unit to heat other lotions, creams or oils.
3. This paraffin wax is reusable. Simply leave the wax in the unit and reuse.
4. Hand, foot or elbow should be lowered carefully into the paraffin, all the way to the wrist or ankle, with the same position being maintained throughout to avoid breaking the paraffin seal.
5. To maintain a constant temperature, always keep the Lid on the unit when not in use.

USER MAINTENANCE INSTRUCTIONS

This unit contains no user serviceable parts. If damaged, refer to the warranty section of this manual for information or return and repair.

CLEANING

NOTE: *Our paraffin wax is reusable, however, if more than one person is using the unit you may wish to replace the wax. Also, it is recommended that you periodically replace the unused wax and clean the unit.*

1. Unplug unit and remove Lid. Allow paraffin wax to solidify overnight.
2. Push down on one side of paraffin wax block to dislodge. **DO NOT USE SHARP OBJECTS TO REMOVE WAX, THIS MAY SCRATCH OR DAMAGE THE UNIT.**
3. If the wax does not easily dislodge, run unit for 10-15 minutes to loosen wax.
4. Discard wax.
5. Wipe the inside of the unit with a soft, slightly damp cloth. Do not pour water into the unit.

STORAGE

Use the Cord Wrap storage area on the bottom of the unit. Carefully wrap cord inside the Cord Wrap Tabs and down through the Cord Slot. Store in its original box or in a cool, dry area. Do not pull on Cord.

LIMITED ONE YEAR WARRANTY

Warranty: This Relaxor® product is warranted by Salton, Inc. to be free from defects in materials or workmanship for a period of (1) year from the original purchase date. This product warranty covers only the original consumer purchaser of the product.

Warranty Coverage: This warranty is void if the product has been damaged by accident in shipment, unreasonable use, misuse, neglect, improper service, commercial use, repairs by unauthorized personnel, normal wear and tear, improper assembly, installation or maintenance abuse or other causes not arising out of defects in materials or workmanship. This warranty is effective only if the product is purchased and operated in the USA, and does not extend to any units which have been used in violation of written instructions furnished with the product or to units which have been altered or modified or to damage to products or parts thereof which have had the serial number removed, altered, defaced or rendered illegible.

Implied Warranties: ANY IMPLIED WARRANTIES WHICH THE PURCHASER MAY HAVE ARE LIMITED IN DURATION TO ONE (1) YEAR FROM THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

Warranty Performance: During the above one-year warranty period, a product with a defect will be either repaired or replaced with a reconditioned comparable model (at our option) when the product is returned to the Repair Center, or the purchase price refunded. The repaired or replacement product will be in warranty for the balance of the one-year warranty period and an additional one-month period. No charge will be made for such repair or replacement.

Service and Repair: Should the appliance malfunction, you should first call toll-free 1 (800) 233-9054 between the hours of 9:00 am and 5:00 pm Central Standard Time and ask for CONSUMER SERVICE stating that you are a consumer with a problem. Please refer to model number PWB8FE when you call.

In-Warranty Service (USA): For an appliance covered under the warranty period, no charge is made for service or postage. Call for return authorization 1 (800) 233-9054 .

Out-of Warranty Service: A flat rate charge by model is made for out-of-warranty service. Include \$ 10.00 (U.S.) for return shipping and handling. We will notify you by mail of the amount of the charge for service and require you to pay in advance for the repair or replacement.

For Products Purchased in the USA, but Used in Canada: You may return the product insured, packaged with sufficient protection, and postage and insurance prepaid to the USA address listed below. Please note that all customs duty / brokerage fees, if any, must be paid by you and we will require you to pay the cost of customs duty / brokerage fees to us in advance of our performing any service.

Risk During Shipment: We cannot assume responsibility for loss or damage during incoming shipment. For your protection, carefully package the product for shipment and insure it with the carrier. Be sure to enclose the following items with your appliance: any accessories related to your problem, your full return address and daytime phone number, a note describing the problem you experienced, a copy of your sales receipt or other proof of purchase to determine warranty status. C.O.D. shipments cannot be accepted.

To return the appliance, ship to:

ATTN: Repair Center
Salton, Inc.
708 South Missouri Street
Macon, MO 63552

To contact us, please write to, call, or email:

Salton, Inc.
P.O. Box 1637
Columbia, MO 65205-1637
1-800-233-9054
E-mail: Salton@Saltonusa.com

Limitation of Remedies: No representative or person is authorized to assume for Salton, Inc. any other liability in connection with the sale of our products. There shall be no claims for defects or failure of performance or product failure under any theory of tort, contract or commercial law including, but not limited to negligence, gross negligence, strict liability, breach of warranty and breach of contract. Repair, replacement or refund shall be the sole remedy of the purchaser under this warranty, and **in no event shall Salton be liable for any incidental or consequential damages, losses or expenses.** Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

Legal Rights: This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

For more information on Salton products: visit our website: www.salton-maxim.com

©2001 Salton, Inc. All rights reserved.

Made in China

P/N 60811

9/01