


Relaxor[®]

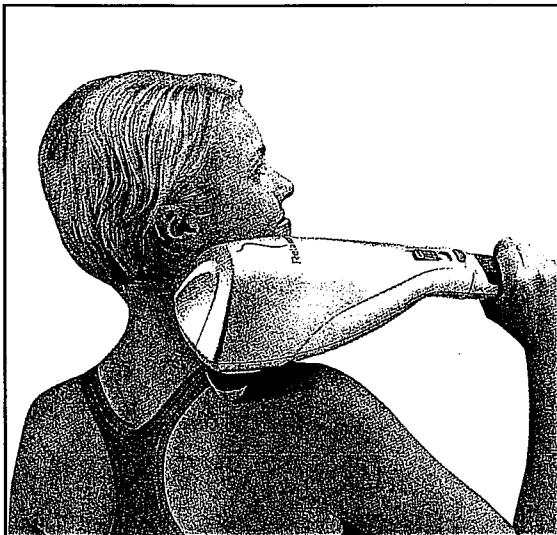

Regenesis[®] *Massager* *with Twin Detachable* *Heads*

**Model
515**

OWNERS MANUAL and OPERATING INSTRUCTIONS

Dear Customer:


Thank you for purchasing the Regenesis[®] Massager with Twin Detachable Heads.

The Regenesis[®] Massager, like our other fine products, is engineered for quality, comfort and safety. It provides penetrating relief and relaxation for sore and tense muscles. This massaging product will also serve as a great addition to your health and fitness program.

In order to enhance your enjoyment and maximize the benefits of your new purchase, please take a few moments to read the following instructions.

Relaxor®

HOW TO OPERATE YOUR REGENESIS® MASSAGER

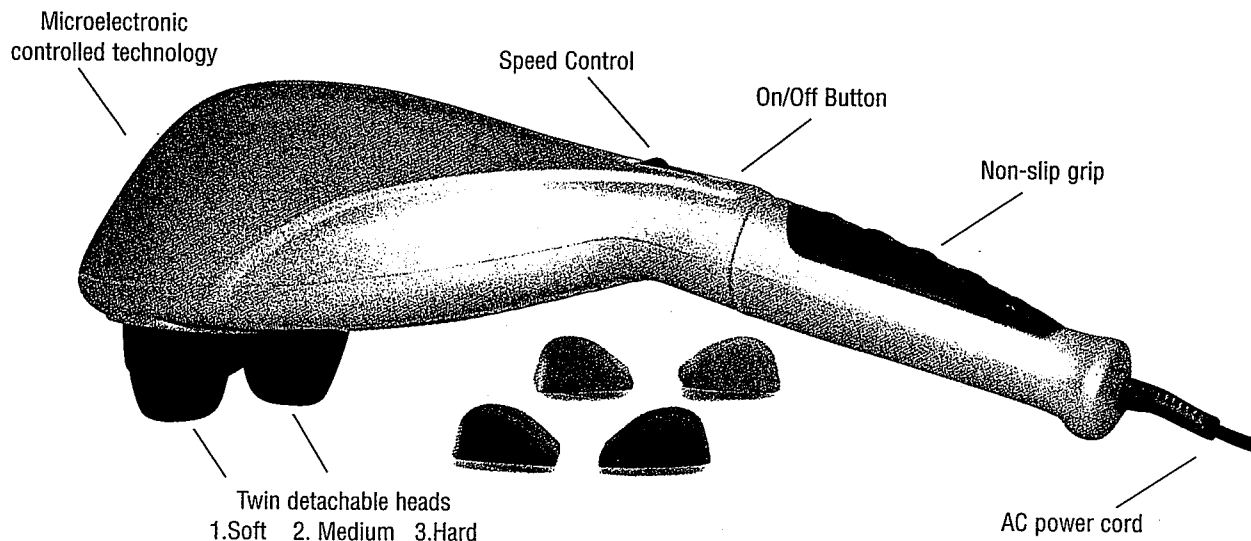

INSTRUCTIONS FOR USE:

1. This product has a polarized plug (*one blade is wider than the other*). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Never use with an extension cord unless plug can be fully inserted. Do not attempt to defeat this safety feature.
2. Always make sure unit is set to "OFF" position before inserting plug into electrical outlet.
3. To turn unit on, press the ON/OFF button, located on back of unit.
4. Adjust the speed control slide switch to the desired position.
5. Use the massager on arms, legs, back or any area you wish to treat. Slowly move the unit around for maximum comfort.
6. When finished, press the power switch to "OFF" and unplug the unit. Store the massager in a dry area.


Changing Massage Attachments:

The Regensis® comes with 3 custom massage attachment sets. Each set is of a different density:

Black = Soft; Medium gray = Medium; Light gray = Hard

Unplug from outlet before putting on or taking off attachments.

To insert, line up cup with straight edge cutaway and push.

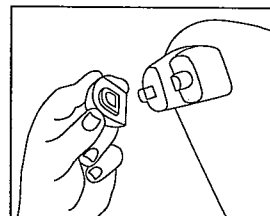

IMPORTANT SAFETY INSTRUCTIONS

Read all instructions before using the Regenesi[®]s Twin Head Massager.

DANGER!

To reduce risk of electric shock:

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not reach for a product that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.
4. Do not place or store product where it can fall or be pulled into a tub or sink.
5. Do not place product in or drop into water or other liquid.

WARNING!

To reduce risk of burns, fire, electric shock, or injury to persons:

1. Your Regenesi[®]s Massager should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. This product has a polarized plug (*one blade is wider than the other*). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Never use with an extension cord unless plug can be fully inserted. Do not attempt to defeat this safety feature.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Do not use near or place in water.
5. Do not place or store where product can fall or be pulled into tub, toilet or sink.
6. Never operate the Regenesi[®]s Massager if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact your dealer for instructions on returning it for examination and repair.
7. Do not use on children or disabled persons.
8. Do not use outdoors.
9. Keep the cord away from heated surfaces.
10. Do not operate where aerosol (*spray*) products are being used or where oxygen is being administered.
11. If you are under a physician's care and/or have a condition requiring the use of any medical device, please consult your health care provider prior to using this product.
12. Do not carry the product by the supply cord or use the cord as a handle.
13. To disconnect, switch all controls to the OFF position, then remove plug from the outlet.
14. Household use only.

OPERATE
YOUR
SYSTEM
ACCORDING
TO THIS
MANUAL

SAVE THESE INSTRUCTIONS!

MAINTENANCE INSTRUCTIONS:

No user maintenance of this product is required. If the System does not function, check to be sure the CORD is plugged into an electrical outlet.


If there is a malfunction of any kind, contact your dealer for instructions on returning it for examination and repair.

Store this product in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects.

Do not wrap the cord around the product. For hand-supported product, untwisting the cord before use.


Relaxor®

OTHER FINE PRODUCTS FROM THE LEADER IN RELAXATION

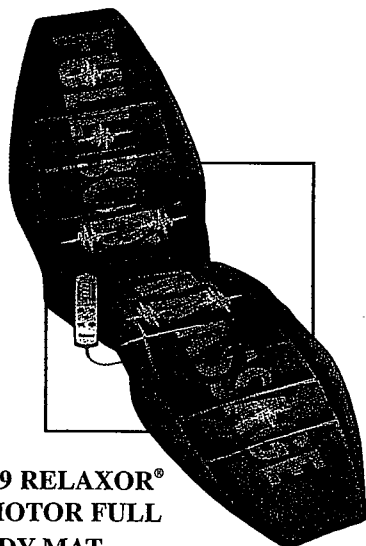


#969 RELAXOR® 7 MOTOR FULL BODY MAT

Enjoy the penetrating caresses of a relaxing massage with soothing heat and the convenience of portability.

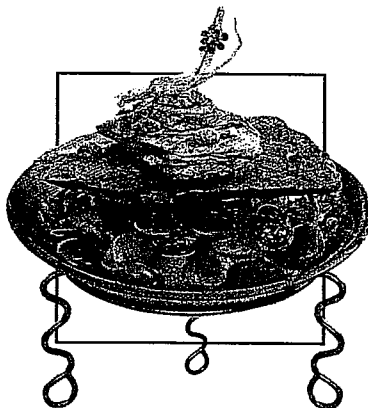

#280 CALMING POOL

Designed to capture the elements of nature by combining water, assorted stones and decorative driftwood to complement any home or office environment. AC powered.

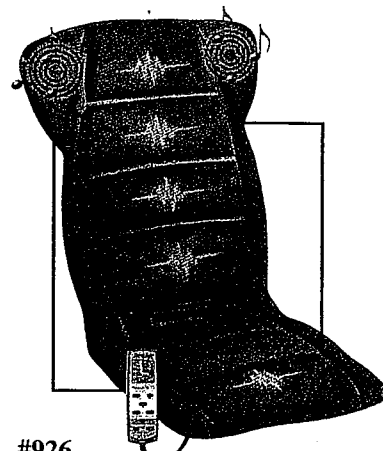


#926 RELAXOR® SEAT TOPPER WITH NATURE SOUNDS

5-Motor Massaging Seat Topper with built-in stereo speakers and Nature Sounds. AC and auto adaptor included for listening in the office, home or auto.


#964 RELAXOR® MASSAGING SEAT TOPPER

12 Motor Massage Seat Topper with Soothing Heat. Features dozens of settings for an all over massage. Great for home, office and auto.

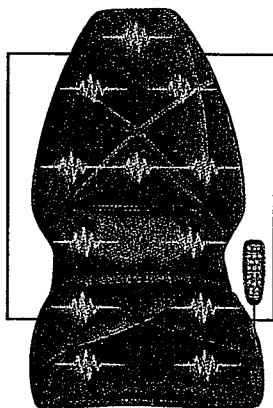

722 STEREO SOUND THERAPY WITH ALARM CLOCK AND AM/FM RADIO

12 soothing environmental sounds to promote rest and relaxation. Full-featured clock and headphone outlet. Customize by mixing any two sounds. Battery and AC powered.

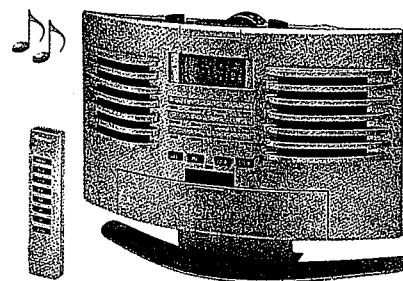

Relaxor®

The Leader in Relaxation™

Relax and Enjoy!

Limited 1-Year Warranty

If your REGENESIS® MASSAGER does not operate properly, or you need assistance for any reason, call:

Customer Service - JB Research, Inc.
Los Angeles Corporate Office • Bellflower, CA

1-800-RELAXOR

Outside U.S. write to:
Customer Service
9062 Rosecrans Ave. Bellflower, CA 90706
© Copyright 1999 JB Research, Inc.

PATENTS PENDING

